

PRODUCT GUIDE

# ECIS


ELECTRIC CELL-SUBSTRATE IMPEDANCE SENSING



Applied  
BioPhysics

[www.biophysics.com](http://www.biophysics.com)

# ECIS Time Course Measurements

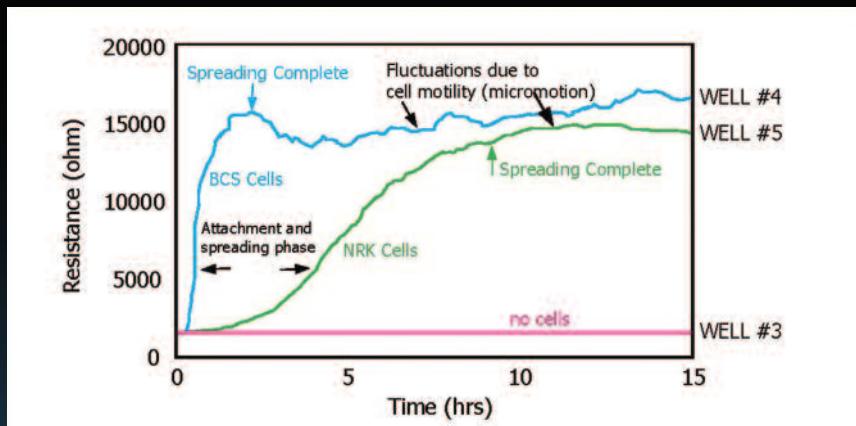

Cells can be sampled as often as several times per sec to as slowly as desired. Each impedance reading is plotted as a point, in Ohms or nanofarads (C), verses time. The total acquisition time is user controlled and can range from a fraction of an hour to several days. The above plot is a measurement of cells attaching and spreading in two different wells of an array. Well #4 are BCS cells and Well #5 are NRK cells. The electrode in each well was sampled at a few second intervals for 15 hours.

Following inoculation at time zero, impedance increases as the cells attach to the electrode and begin spreading. The impedance continues to increase until the cells reach confluence at 2 hours for the BCS cells and 10 hours for the NRK cells. The small fluctuations in the curves are due to micromotion from the constant movement of the monolayer of cells on the electrode.



Applied  
BioPhysics

## TABLE OF CONTENTS

<b>1. What is ECIS</b> .....	2–5
How to Quantify Cell Behavior.....	3
How Frequencies Reveal Cell Behavior .....	4
How Electrode Designs Reveal Aspects of Cell Behavior.....	5
<b>2. Applications</b> .....	6–15
Attachment and Spreading .....	7
Cell Proliferation .....	8
Differentiation and Stem Cell Biology .....	8
Barrier Function .....	9
Signal Transduction .....	9
Cell Invasion .....	10
Cell Toxicity .....	11
Automated Cell Migration .....	12
In SITU Electroporation & Monitoring .....	13
Inflammation .....	14
Flow Based Assays.....	15
<b>3. Instruments and Accessories</b> .....	16–19
ECIS ZØ .....	17
16 and 96 Well Array Stations.....	17
TEER24.....	18
CP96 Cell Proliferation.....	18
Array Stage Holder.....	19
P Flow Peristaltic Pump.....	19
CO2 Tissue Culture Incubator.....	19
<b>4. Arrays</b> .....	20–27
96W Standard Arrays .....	22
8W Standard Arrays .....	23
Specialty Arrays.....	24–25
Flow Arrays .....	26–27
Institutions That Use ECIS .....	28–29
International Distributors.....	30
Request a Demo.....	32

# 1 what is ECIS


ECIS® (Electric Cell-substrate Impedance Sensing) is a real-time, impedance-based method to study many of the activities of cells when grown in tissue culture. These include morphological changes, cell locomotion, and other behaviors directed by the cell's cytoskeleton. Impedance-based cell monitoring technology was invented by Drs. Ivar Giaever and Charles R. Keese who formed Applied BioPhysics, Inc. to commercialize ECIS® and other biophysical technologies.

The ECIS® approach has been applied to numerous investigations including measurements of the invasive nature of cancer cells, the barrier function of endothelial cells, in vitro toxicity testing as an alternative to animal testing, and signal transduction involving GPCR's for modern drug discovery.



Applied  
BioPhysics


## HOW TO QUANTIFY CELL BIOLOGY

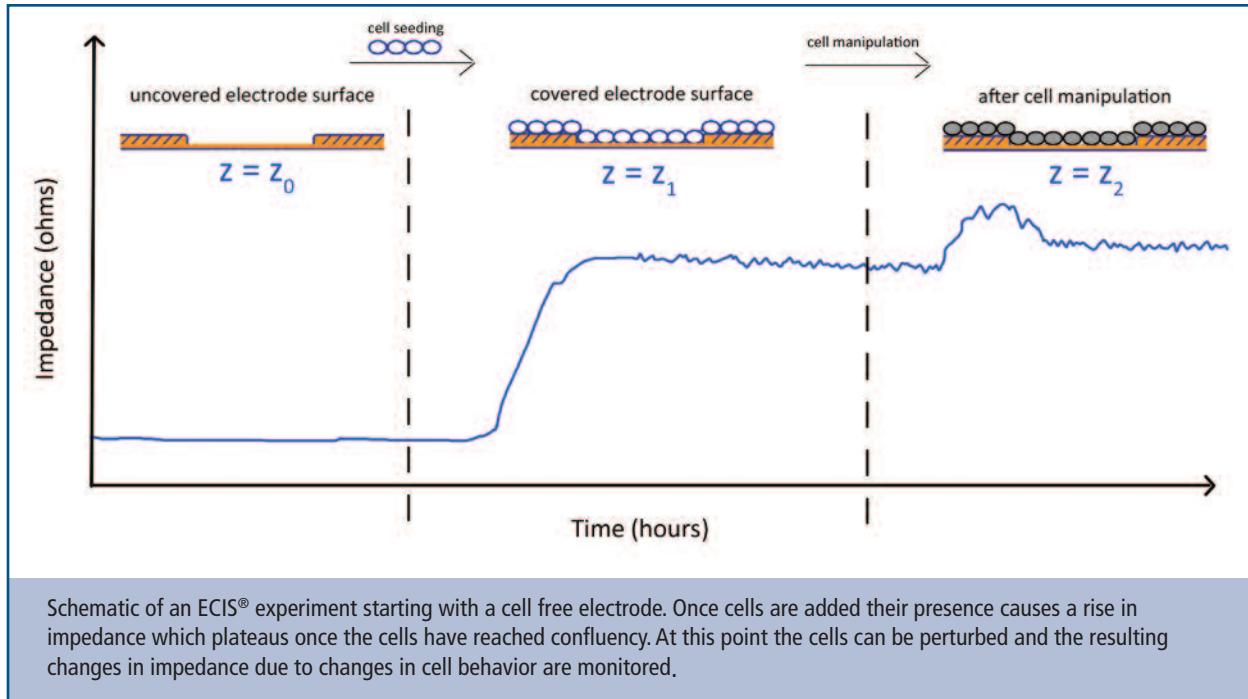

Cell function modulates cell morphology. ECIS® is capable of detecting and quantifying morphology changes in the sub-nanometer to micrometer range. In ECIS® a small alternating current ( $I$ ) is applied across the electrode pattern at the bottom of the ECIS® arrays (direct current cannot be used). This results in a potential ( $V$ ) across the electrodes which is measured by the ECIS® instrument.

The impedance ( $Z$ ) is determined by Ohm's law  $Z = V/I$ . When cells are added to the ECIS® Arrays and attach to the electrodes, they act as insulators increasing the impedance. As cells grow and cover the electrodes, the current is impeded in a manner related to the number of cells covering the electrode, the morphology of the cells and the nature of the cell attachment.


When cells are stimulated to change their function, the accompanying changes in cell morphology alter the impedance. The data generated is impedance versus time.

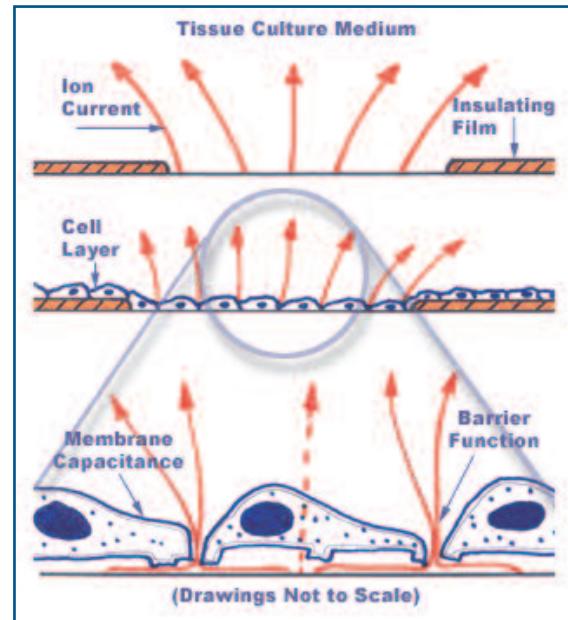


ECIS® electrodes are defined by small openings in an insulating film covering the gold surface. Once cells fill these openings the current flow from the gold surface is restricted.

## HOW FREQUENCIES REVEAL CELL BEHAVIOR

To understand why AC frequency is important in ECIS® we have to consider how frequency affects the current paths of cell-covered electrodes. (Note: the total current is maintained constant and voltage changes are measured.) At relatively low frequencies (< 2,000Hz) most of the current flows in the solution channels under and between adjacent cells (red lines).

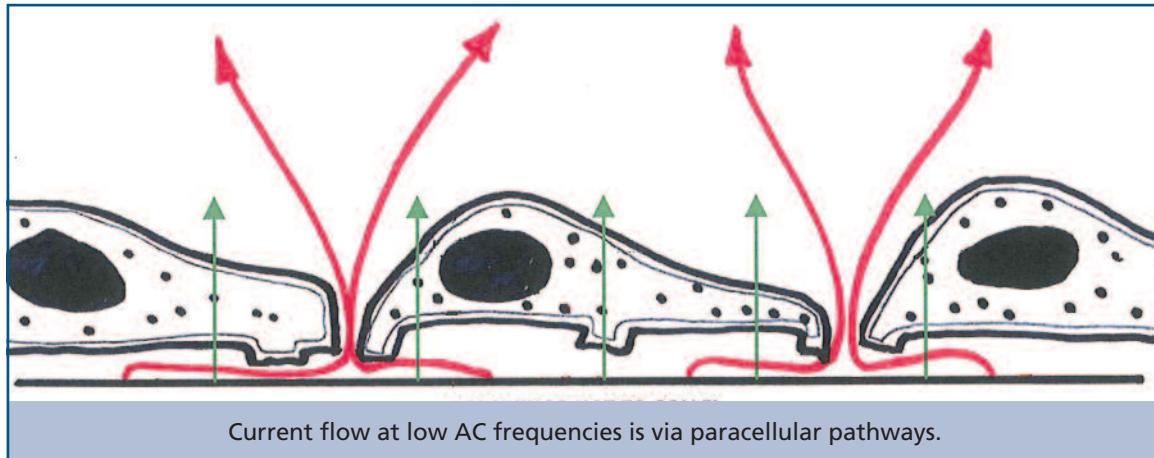

At higher frequencies (> 40,000 Hz) more current now capacitively couples directly through the insulating cell membranes (green lines).

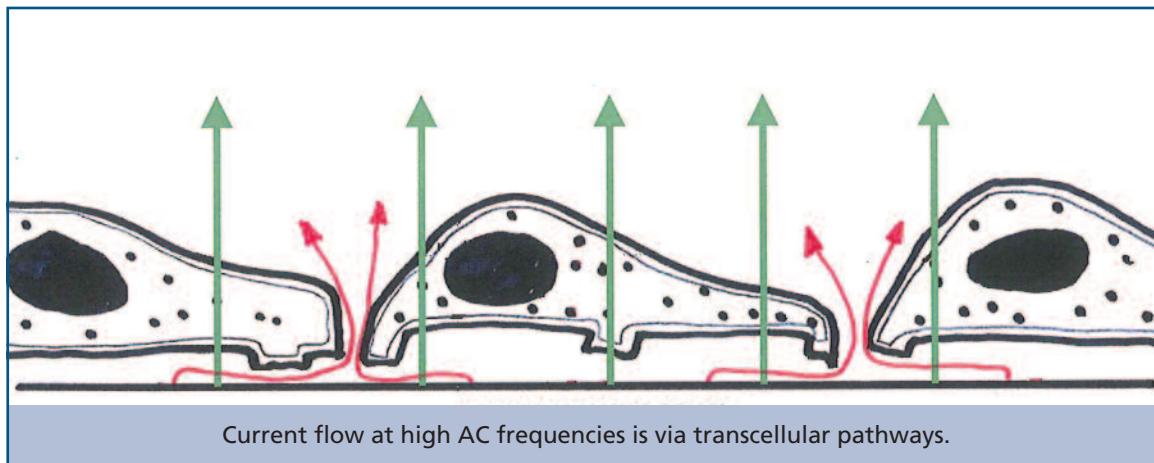


The high frequency impedance is more affected by cell-coverage, whereas the low frequency responds more strongly to changes in the spaces under and between the cells.

With the more advanced ZΘ instrument, where the impedance is broken down into its components (resistance and capacitance), quantitative information about the cells can be obtained by modeling (Giaever and Keesee PNAS 1991).

Using impedance data at multiple AC frequencies the ECIS® model calculates time course changes in:

- The barrier function (permeability) of the cell layer
- The degree of constricted flow of current under the cells
- The cell membrane capacitance

## HOW ELECTRODE DESIGNS REVEAL ASPECTS OF CELL BEHAVIOR

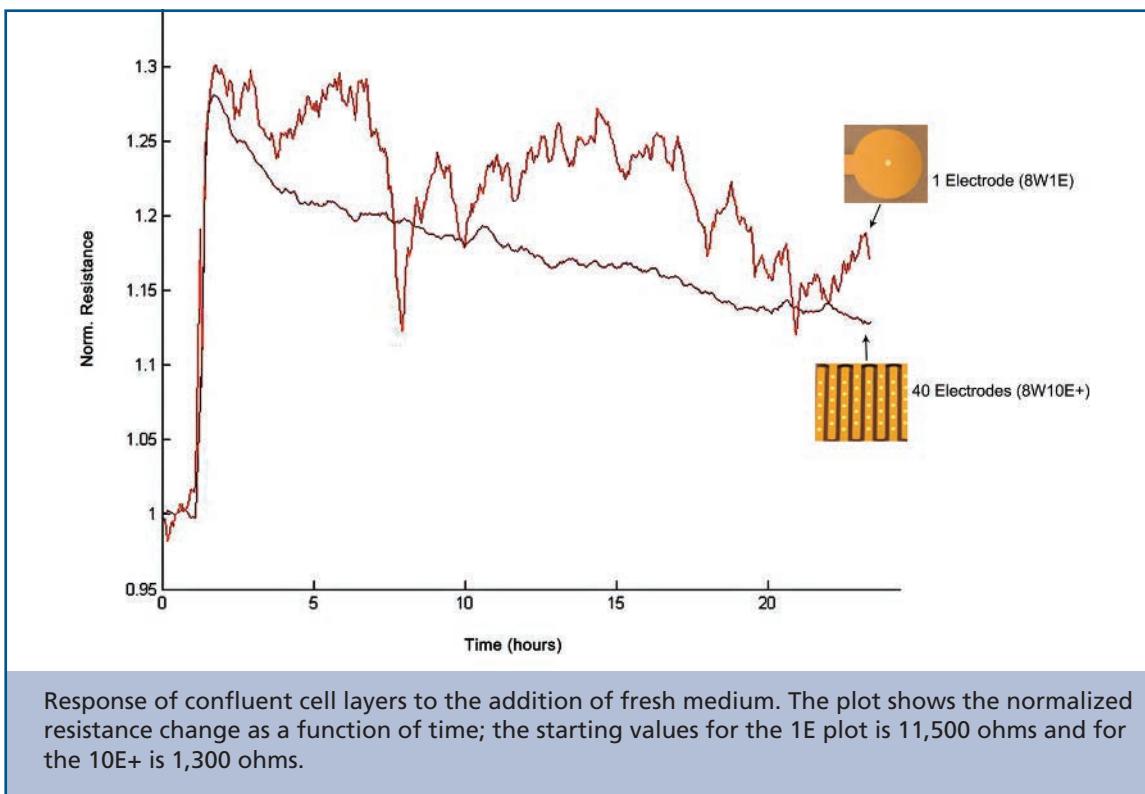

### Small Electrodes

Small electrodes (1E, 10E, 10E+ type arrays) and their layout within the wells ensure that all current passes through the cell monolayer. This allows the ability to analyze data with the ECIS® modeling software to determine barrier function, cell membrane capacitance as well as the spacing between the cell basal membrane and electrode.

Keeping the total surface area of the electrodes small also allows for a relatively low AC current to generate the large electric field necessary to either electroporate or kill the cells in migration experiments.

Small electrodes also provide the ability to monitor the uncorrelated nano-scale morphological changes of individual or small populations of cells (<100), while larger or multiple electrodes provide the averaged morphological response of many cells (1000+).

### Large Electrodes

Some experimental protocols, such as cell proliferation, require sparse inoculations leading to a variance of cell density at the bottom of the well. Large electrodes (CP Array) or a large collection of small electrodes (10E+ Array) increases the sampling size resulting in less variability.


# 2 applications

By quantifying cell behavior ECIS® characterizes the life cycle of the cell: attachment and spreading of cells, growth into a confluent monolayer, the dynamics of the monolayer, reactions to stimuli, and finally cell death.



Applied  
BioPhysics



## ATTACHMENT AND SPREADING

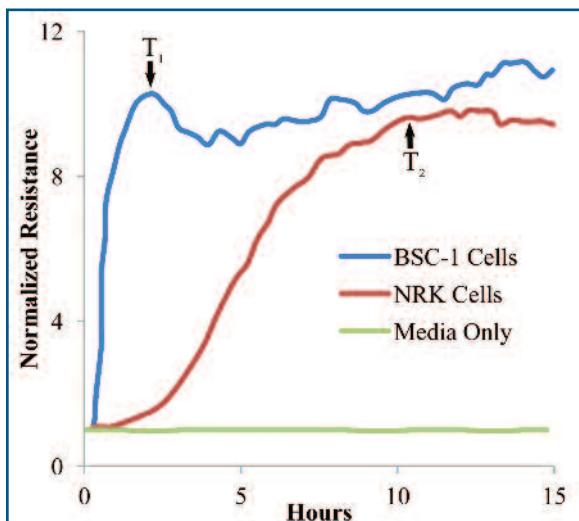

ECIS® graph of the attachment and spreading of BSC-1 and NRK cells onto ECIS 1E type electrodes. Cells were seeded at density of 105 cells/cm<sup>2</sup>, and a confluent monolayer was achieved at T<sub>1</sub> and T<sub>2</sub> for BSC-1 and NRK cells respectively. The fluctuations that exist after these time points are due to nano-scale cellular activity.

ECIS® experiments start with freshly seeded cells attaching to the ECIS® electrodes at the bottom of the wells. As cells attach to the surface they begin to spread, increasing the amount of cell area in contact with the ECIS® electrode and the measured impedance. ECIS® provides a continuous real-time measurement of this process, quantifying both the rate of change and the final value of the impedance. Published examples of attachment assays include dependencies upon ECM proteins, genetic manipulation of junction and signal transduction proteins, and binding competition with anti-integrin antibodies or the tetrapeptide RGDS.


Traditional "counting attached cells assays" can only quantify the number of cells attached to any ECM coating. ECIS® assays give feedback on the strength of the attachment of the cells to the ECM. The transparent nature of the electrode mean results can be normalized to resistance per cell or per cell area by means of optical inspection. Longer term assays can probe not just the effect of ECM proteins on attachment, but functional consequences of ECM derived from different sources (see B below).

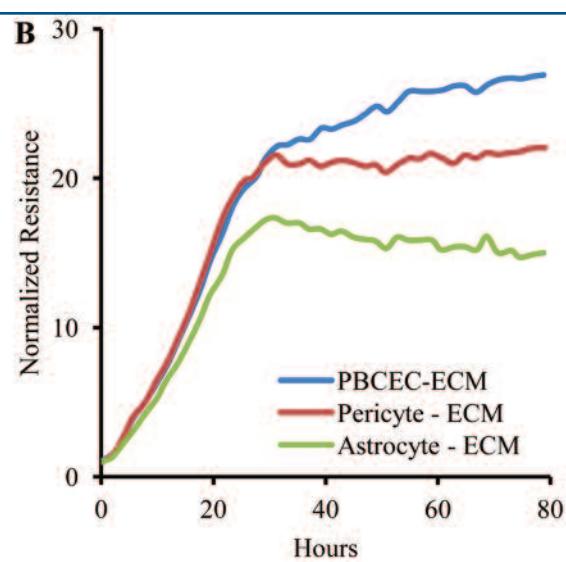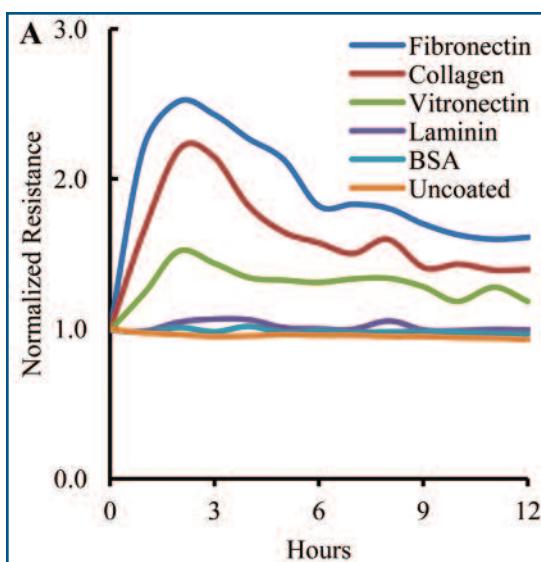

Graphs depicting the attachment and spreading dynamics of (A) renal vascular smooth muscle cells and (B) porcine brain capillary endothelial cells. In (A) ECIS® arrays were pretreated with the indicated matrix proteins prior to cell seeding. In (B) prior to recording, pericytes, astrocytes, or cerebral endothelial cells were allowed to grow to confluence in ECIS® arrays. The cells were then removed by trypsinization and porcine brain endothelial cells introduced to the wells and their attachment and spreading dynamics were measured. Data of (A) derived from Balasubramanian, L. et al., 2008 Am. J. Physiol. Cell Ph. 295:C954 and (B) from Hartmann, C. et al., 2007 Exp. Cell Res. 313:1318.

## CELL PROLIFERATION

As cells proliferate two factors act to change the impedance: cell number and cell morphology. In most instances the cells grow asynchronously and the impedance gradually increases until a maximum when cells become confluent. The impedance change is approximately linear with cell number while the cells are sub-confluent. If the cells are synchronized, then the progression of cell morphological changes associated with the cell cycle can be visualized as a series of peaks and valleys in the impedance plot. Compounds affecting cell growth can be introduced before or after cells have attached to distinguish changes in growth from the ability of cells to attach to the substrate.

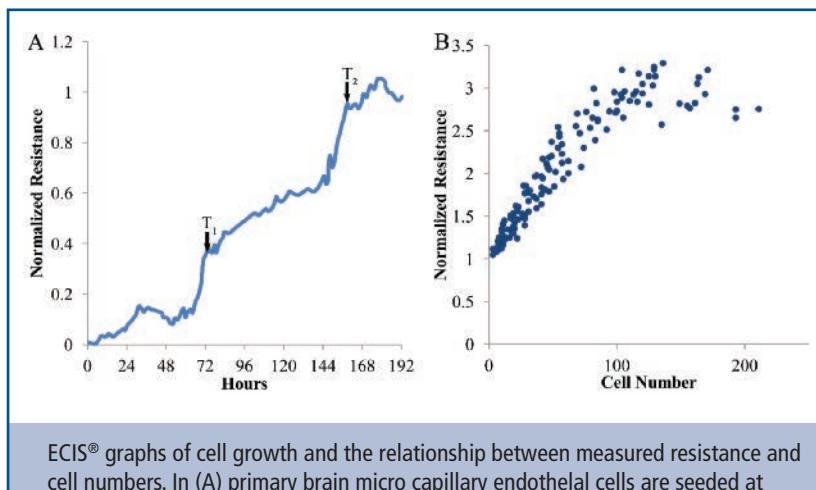


ECIS® graphs of cell growth and the relationship between measured resistance and cell numbers. In (A) primary brain micro capillary endothelial cells are seeded at time zero and continuously measured over the ensuing week. Plateaus indicated at T<sub>1</sub> and T<sub>2</sub> represent time points at which the cell population has doubled. (B) A linear correlation of resistance with the number of cells on the electrode exists below a saturation density of approximately 100 cells per electrode. Data derived Bernas, M.J. et al., 2010 Nat. Protocols 5:1265 (A) and from Xiao, C. & Luong, J.H.T., 2003 Biotechnol. Progr. 19:1000

## DIFFERENTIATION AND STEM CELL BIOLOGY


When cells differentiate they change their behavior allowing ECIS® to follow the events of cell differentiation. While most tools available to characterize stem cells preclude their further use, the label-free non-invasive nature of ECIS® allows for subsequent use of characterized stem cells. Thus stem cells can be selected for functional abilities before their use in research or therapy.

This ability to characterize cells by their function allows ECIS® to be a powerful tool in phenotypic drug discovery. It use is not limited to changes in phenotypes but in also assuring the quality of cell stocks or assurance that differentiated cells have not reverted

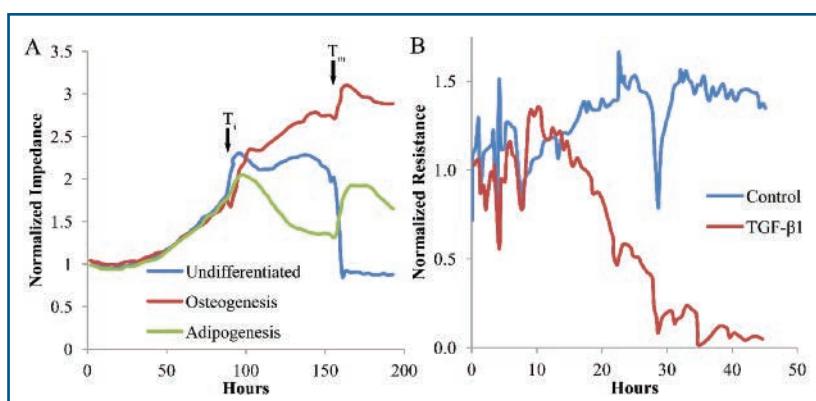

ECIS® graph of differentiating cells. In (A) adipose-derived stem cells are monitored before and after differentiation is induced at T<sub>1</sub>. The distinct differences between the differentiated cell types persists beyond the first media change at T<sub>m</sub>. In (B) NMuMG (murine mammary gland) cells are monitored after the addition of 5ng/ml TGF-β1 (red line) at 0 hours. The decline in normalized resistance between 10 and 30 hours correlates with the transition of the NMuMG cells from having endothelial characteristics to mesenchymal. Data derived from Bagnaninchi, P.O. & Drummond, N., 2011 PNAS 108:6462 (A) and Schneider, D., et al., 2011 BBA 1813:2099 (B).


## BARRIER FUNCTION

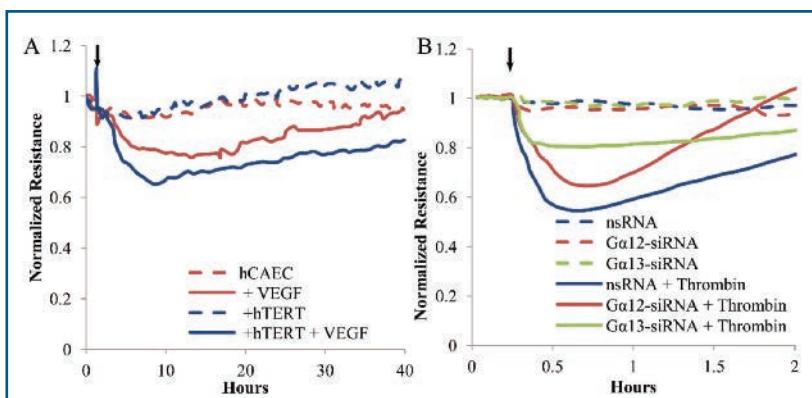

ECIS® graphs of barrier function in HCAEC (A) and human pulmonary aortic endothelial cells (B). In (A) over-expression of human telomerase reverse transcriptase enhances the ability of VEGF to reduce barrier function in human coronary artery endothelial cells. In (B) HPAEC were incubated with siRNA to Gα12 and Gα13 or treated with nonspecific RNA and were stimulated with thrombin (50 nM) at the time marked by the arrow. Data derived from Baumer, Y. et al., 2011 Exp. Bio. Med. 236:692 (A) and (B) Birukova, A. et al., 2004 FASEB J. 18:1879

Epithelial cells and endothelial cells regulate the passage of molecules across cell layers. Diseases, especially vascular disease, occur when this function is impaired. Passage of molecules across an endothelial or epithelial cell layer occurs in two ways; actively by transport through the cell or passively by diffusion in the para-cellular space. ECIS® measurements of the resistive portion of impedance at frequencies below 5kHz are very sensitive to changes in the barrier function of these cell types. ECIS® has been used to demonstrate the effects of many regulating molecules including VEGF, thrombin, TNFalpha, and sphingosine-1-phosphate.

## SIGNAL TRANSDUCTION

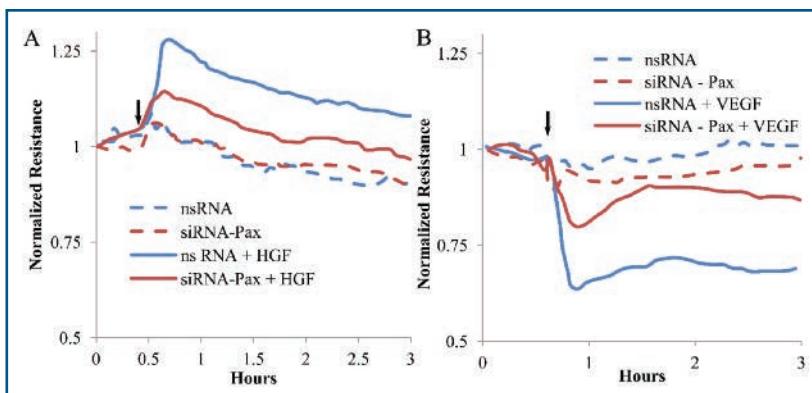


ECIS® graphs for studying the intersecting signal transduction pathways of hepatocyte growth factor (HGF) and vascular endothelial growth factor (VEGF) in pulmonary endothelial cells. (A) The barrier enhancing effect of HGF (added at arrow) is attenuated by the loss of paxillin due to siRNA knockdown. (B) Barrier loss due to VEGF addition (arrow) is also attenuated by the loss of paxillin due to siRNA. Data derived from Birukova, A. et al., 2009 Am. J. Resp. Cell Mol. 40:99


Many compounds of interest are those molecules that specifically bind to cell-surface receptors. When ligands bind to their receptors they initiate signal transduction pathways resulting in a cellular response. ECIS® monitors cellular responses to ligands directly, avoiding false positives, providing reliable information regarding the efficacy of compounds. ECIS® is especially useful to monitor the signal transduction pathways activated by G protein coupled receptors (GPCR). GPCR activation, regardless of the second messenger, results in alterations of the cell's cytoskeletal elements, causing morphological changes. This is precisely the type of cell behavior detected in real time and with great sensitivity by ECIS®.

## CELL INVASION

By quantifying cell behavior, ECIS® can give new insight into how invasive cells and pathogens cross endothelial and epithelial monolayers. By simultaneously monitoring both barrier function and cell viability, ECIS® can distinguish between transmigration mechanisms that leave the monolayer intact from those that disrupt the cell layer. Published examples include metastatic cell and leukocyte trans-endothelial migration, as well as the migration of pathogens such as yeast, anthrax, streptococcus, plasmodium, trypanosomes, and spirochetes.

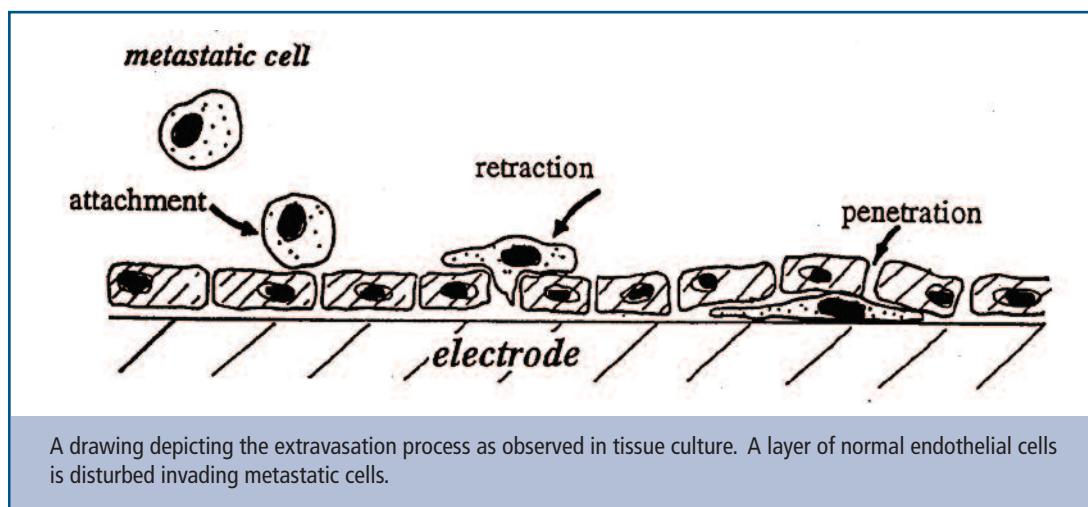


A drawing depicting the extravasation process as observed in tissue culture. A layer of normal endothelial cells is disturbed by invading metastatic cells.

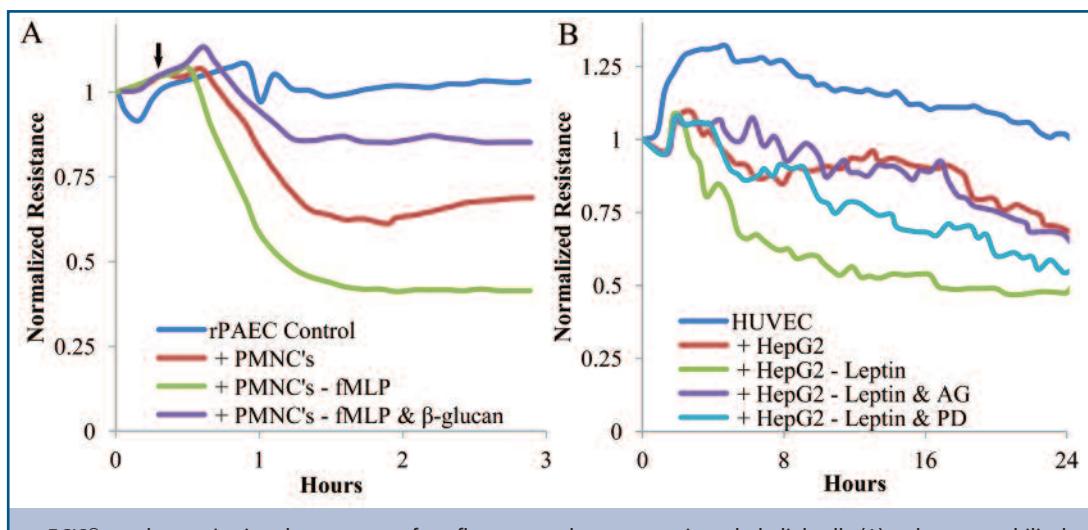

ECIS® graphs monitoring the response of confluent rat pulmonary aortic endothelial cells (A) or human umbilical vein endothelial cells (B) after being challenged by non-adherent cells. In (A) peripheral mononuclear cells are untreated, treated with the chemoattractant fMLP, or treated with fMLP and  $\beta$ -glucan and then added to the confluent rPAEC layer (arrow). In (B) HepG2 hepatocarcinoma cells are untreated, treated with leptin or cotreated with leptin and the JAK/STAT inhibitor AG490 (Calbiochem) or the MAPK inhibitor PD098059 (Sigma). After treatment, the carcinoma cells are added onto the confluent HUVEC layer just prior to recording. Data from (A) Tsikitis, V.L., et al., 2004 J. Immun. 173:1284 and (B) Saxena, N.K. et al., 2007, Cancer Res. 67:2497.


## CELL TOXICITY

Traditionally, toxicity assessments involve animal studies that are both time-intensive and costly. The ECIS® system has been used specifically to assess the cytotoxicity of a variety of toxicants. ECIS-based toxicity tests are far superior to simple cell death assays, because cell function is also monitored. This is important as a change in cell function is not necessarily a cytotoxic event but could be toxic to the individual organism. Cells can be treated in suspension with a toxicant, and then their ability to adhere and grow on the ECIS® electrodes can be monitored. Alternatively, the impedance of established cell monolayers can be monitored after toxicant addition. Analysis of cell micromotion has been shown to enhance the detection of some toxicants by a factor of 10.

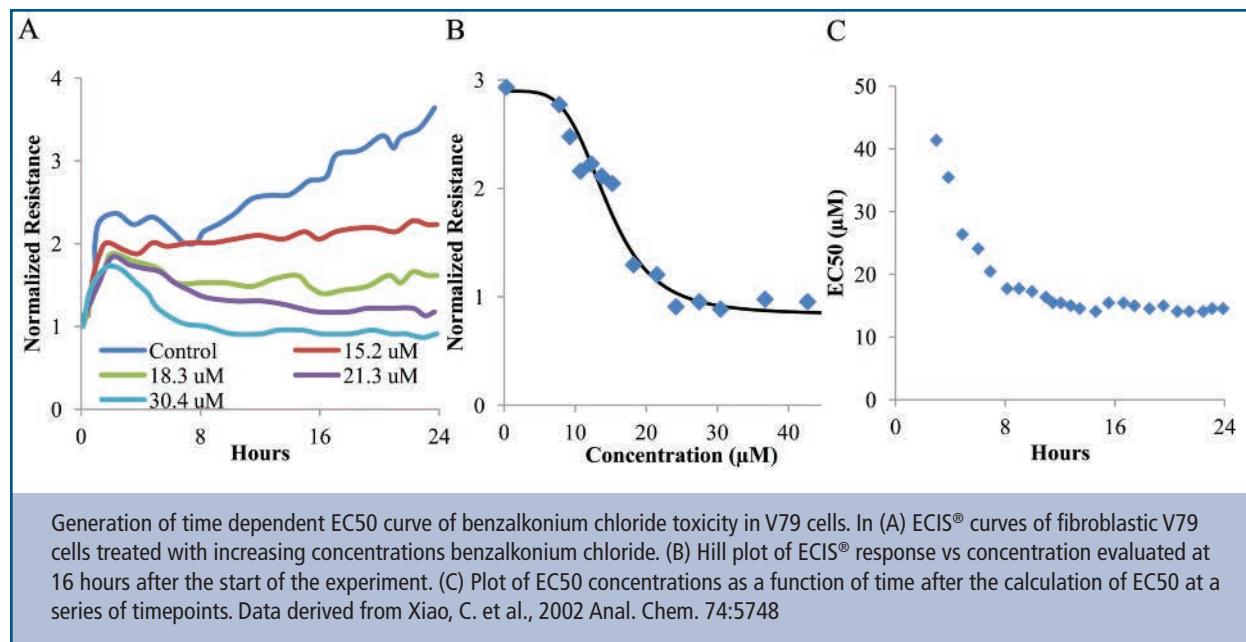


## AUTOMATED CELL MIGRATION

ECIS instruments include an elevated field mode allowing for electroporation and wounding. The ECIS® wound is precisely defined, as it includes only those cells on the electrode. Additionally, with ECIS® the ECM protein coating is not scraped off and is unaffected by the current.

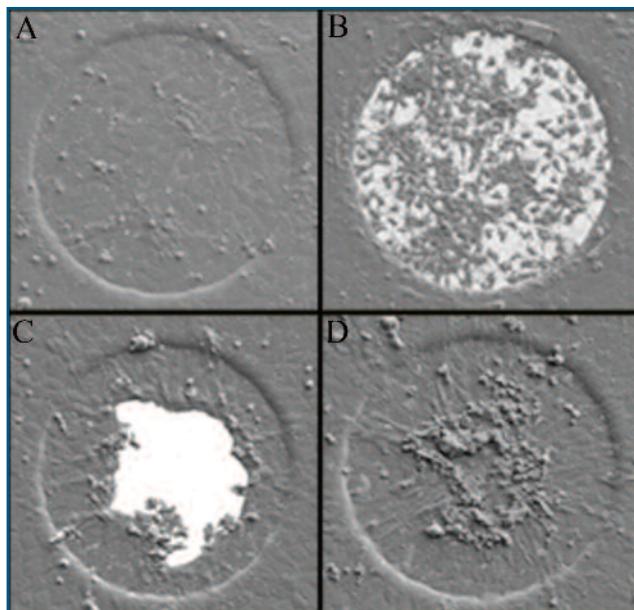


Scanning electron micrographs of the ECIS electrode at time points just prior to (A), just after (B), 4 hours after (C), and 8 hours after (D) the application of a high field pulse across the ECIS electrodes.

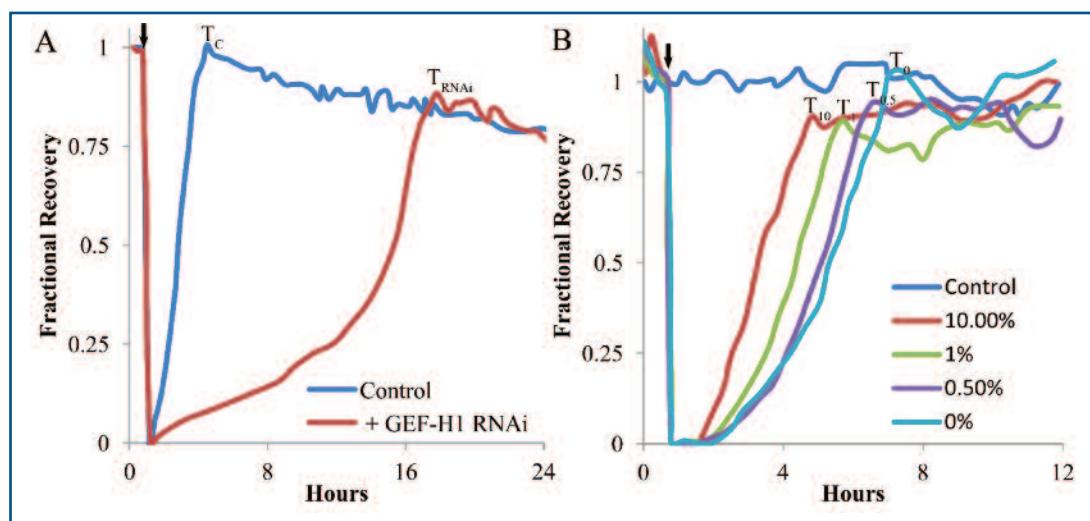


ECIS® graphs of recovery after wounding by applying a high field pulse across the ECIS® electrodes (arrows). In (A) retinal pigment epithelial cells lacking the guanine nucleotide exchange factor GEF-H1 show a prolonged recovery time (T<sub>RNAi</sub>) as compared to their non-treated controls (T<sub>C</sub>). In (B) BSC-1 cells were grown to confluence, and 20 hours prior to the experiment the complete media was replaced with media containing the indicated amount of serum. Recovery times (T<sub>10</sub>, T<sub>1</sub>, T<sub>0.5</sub>, T<sub>0</sub>) were inversely correlated with amount of serum. Control cells were maintained in 10% serum and did not receive a high field pulse. Data derived from (A) Tsapara, A. et al., 2010 Mol. Bio. Cell 21:860 and (B) Keese, C.R., et al., 2004 PNAS 101:1554.

## IN SITU ELECTROPORATION & MONITORING

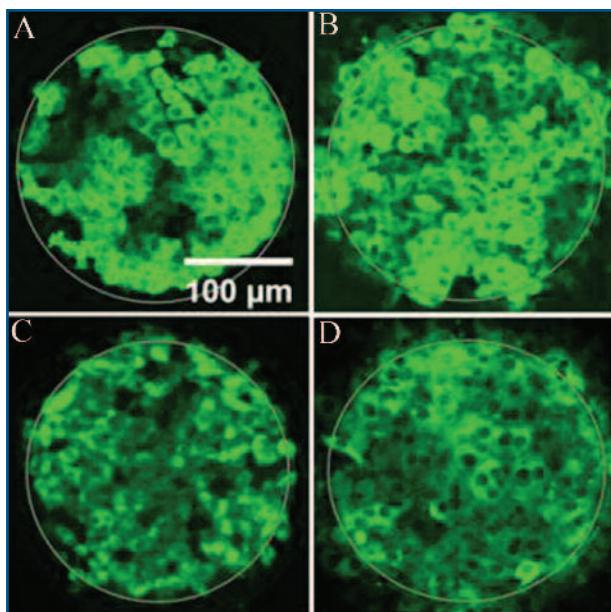

Direct fluorescent imaging of the ECIS® electrode after 250kD FITC-labeled Dextrans were introduced into (A) NRK, (B) HEK-293, (C) CHO, and (D) NIH-3T3 cells by ECIS® in situ electroporation.


ECIS® instruments include an elevated field mode allowing for electroporation and wounding. ECIS® in situ electroporation of adherent cells is more convenient and less disturbing for anchorage-dependent cells than standard electroporation techniques. In addition, it offers the added benefit that the cellular response can be studied immediately after pulsing. Membrane-impermeable toxins, inhibitors or other bioactive compounds can be introduced into the cytoplasm of adherent cells. Monitoring the cell response in real-time, from minutes to weeks, opens up new strategies to study such complex processes as signal transduction, cell differentiation, cell proliferation and apoptosis.

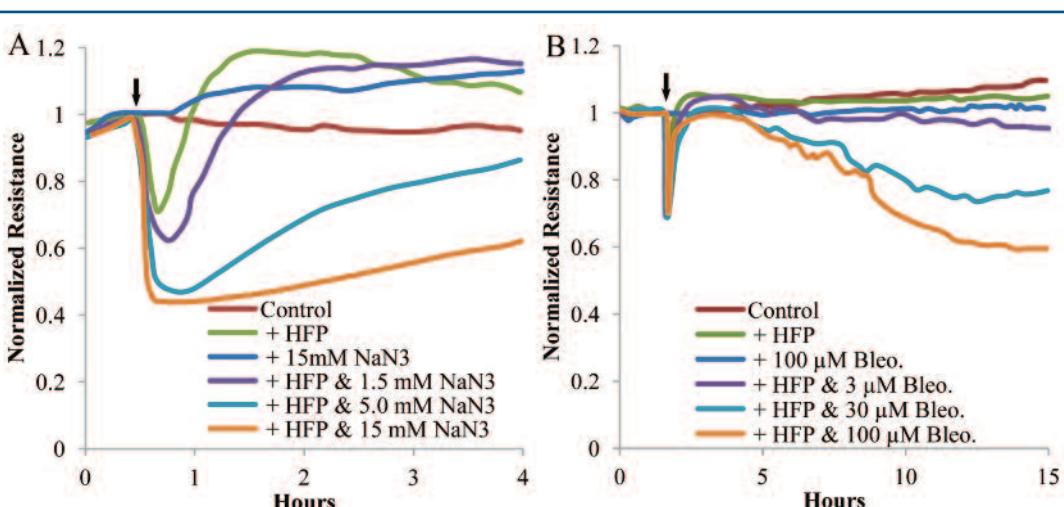

NRK cells incubated with different concentrations of sodium azide (A) or bleomycin (B) at the concentrations indicated and subject to in situ electroporation (Arrowhead) by a high field pulse (HFP). Control cells remained entirely untreated, whereas +HFP cells were electroporated only. A third control set was incubated with the highest concentration of toxin but never electroporated. Derived from Stolwijk, J.A. et al., 2011 Biosens. Bioelectron. 26:4720

## INFLAMMATION

Inflammation is the body's response to pathogens, adverse stimuli such as toxins or ischemia, and physical injury. During an inflammatory response cytokines and interleukins are secreted to guide immunological cells to the site of infection or wounding. Key targets of these molecules are epithelial and endothelial cells which are often the site of injury or need to allow the passage of immune cells to the site of injury. Acute inflammation is generally healing in nature and generally understood. Chronic inflammation is the lack of proper immune regulation and is the core dysfunction for many diseases, including asthma, arthritis, inflammatory bowel disease, cancer, and allergies. Chronic inflammation and its causes are less well understood. ECIS offers a number of cell based assays used to study the inflammatory process. ECIS recovery-after-wounding assays allow for the discovery of molecules which aid in the process of tissue repair. ECIS barrier function assays specifically measure the response of epithelial and endothelial cells to secreted cytokines and can give indirect information about the binding of immune cells to the epithelium or endothelium. Also our newly introduced trans-Filter adapters can be used to study the 3D migration of cells across a Matrigel® layer.


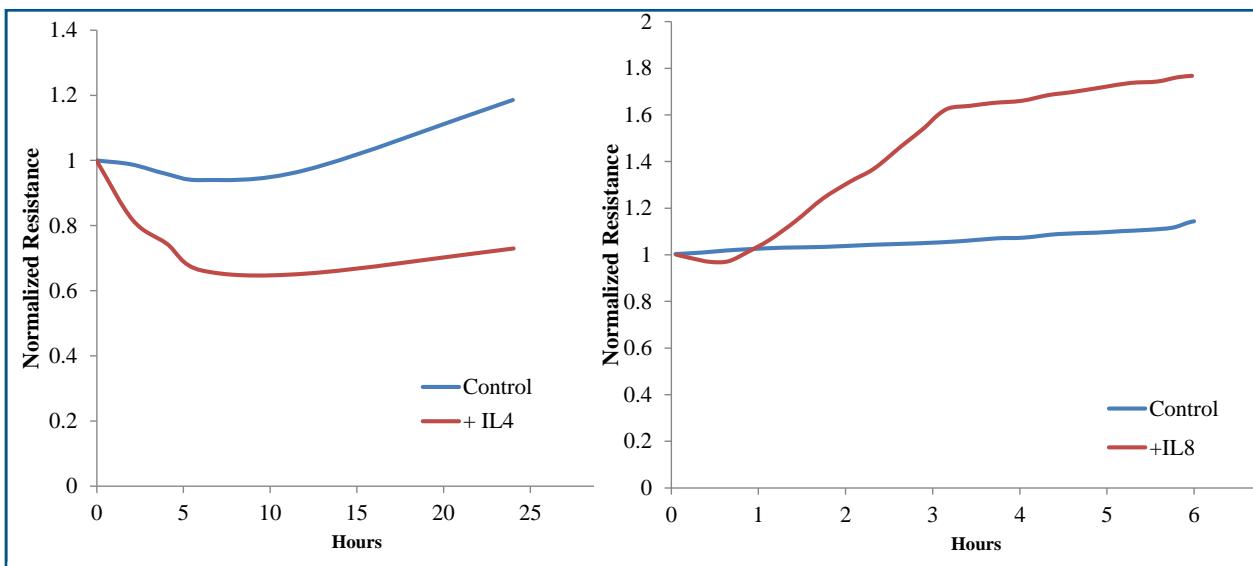


ECIS assays for immunology. In an ECIS barrier function assay (A) bronchial airway epithelial cells respond to the Th2 cytokine il-4 by a drop in resistance. In an ECIS recovery-after-wounding assay(B) keratinocytes demonstrate a faster rate of recovery when treated with the cytokine il-8. Data from (A) derived from Ramirez-Icaza, G. et al. 2004 J. Clin. Immunol. 24:426 and (B) from Jiang, W. G., et al. 2012 Exp. Ther. Med. 3:231.

## FLOW BASED ASSAYS

In their natural environment, endothelial cells are constantly exposed to physical and biochemical stimuli that can alter cell permeability. Laminar shear stress due to blood flow is a principal regulator of systemic endothelial cell gene expression, morphology, and the production of soluble mediators. Its importance is highlighted by pathological processes associated with reduced or absent laminar shear stress, including atherosclerosis. Endothelial transport of macromolecules has been shown to be responsive to flow shear stress, hydrostatic pressure, thermal shock, and agonists such as histamine and thrombin. The ECIS pFlow solutions combines different pump technologies with ECIS flow arrays to allow researchers to study endothelial permeability *in vitro* under complex shear flow conditions.


The pFlow systems use peristaltic pumps to create shear stress conditions within ECIS flow arrays. The ECIS system then monitors continuously the TEER of cell monolayers exposed to the shear stress conditions. This allows for the dynamic changes in TEER to be recorded due to changes in flow rates, addition of vasoactive compounds under flow conditions, or the introduction of secondary cells. Stock ECIS flow arrays create laminar shear stress across the cells. Turbulent flow can be created by customizing the ECIS flow arrays at reasonable costs.

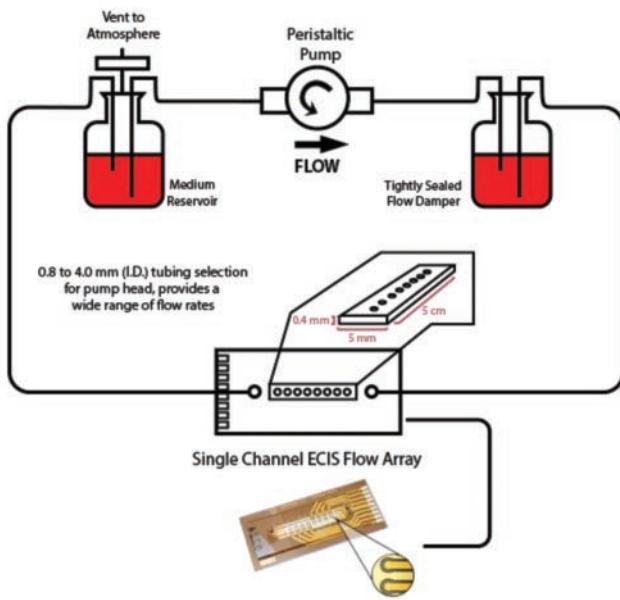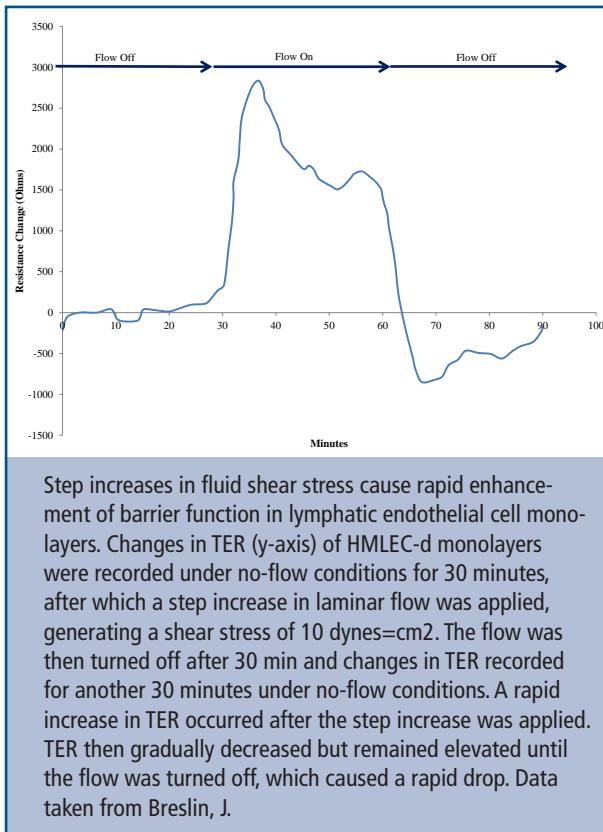


# 3 instruments and accessories

From our all-in-one system to our more function-specific systems, the ECIS line of instruments include different models that range in performance and capability.

Applied BioPhysics has developed a number of accessories that integrate with the ECIS® line of instruments. We offer the following items either directly from Applied BioPhysics or in conjunction with third party vendors.

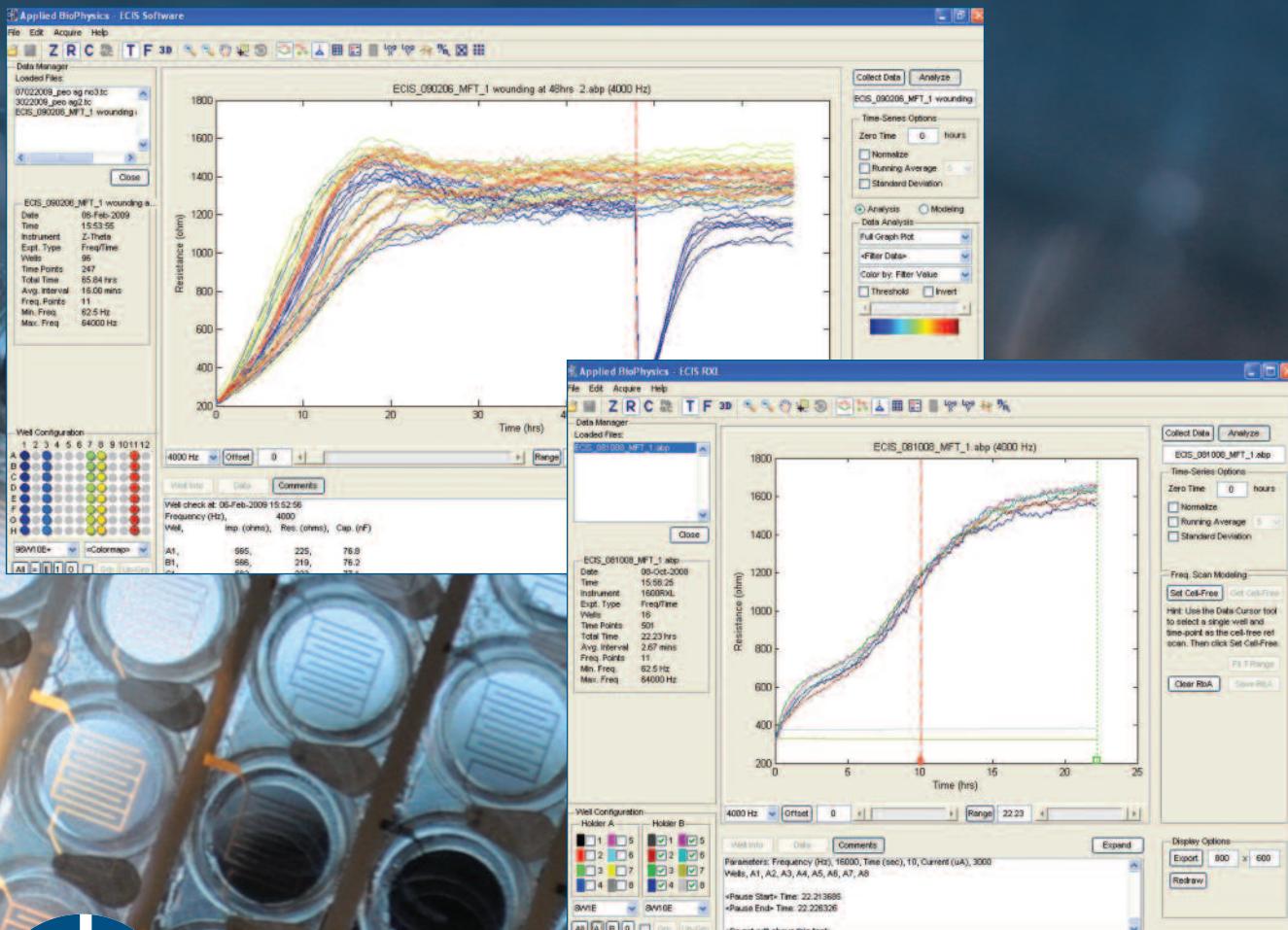

Applied  
BioPhysics

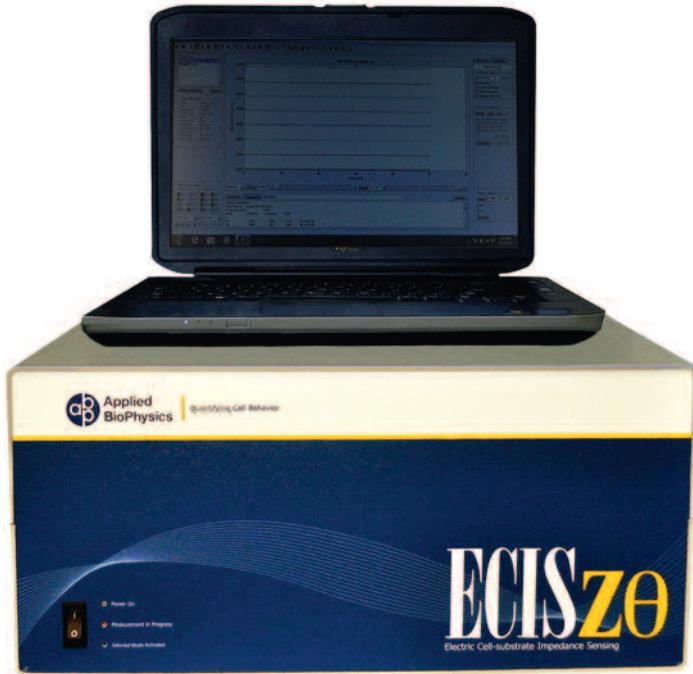

## ECIS® Zθ

The ECIS® Zθ (Z Theta) is our most advanced instrument capable of noninvasively monitoring cell behavior in real-time. The turn-key package consists of a Zθ System Controller, 16 and/or 96 well Array Station, and computer with integrated software running on Windows. The system also includes our exclusive elevated field mode (EFM) to carry out automated cell migration and electroporation studies.

The ECIS® Zθ interprets complex impedance as series resistance and capacitance and reports these values at any AC frequency. Data gathered from confluent cell monolayers at multiple AC frequencies can also be mathematically modeled to present time course changes in three parameters related to actual cell morphology. These parameters include changes in the barrier function (permeability) of the confluent layer, averaged apical and basal membrane capacitance and the close contacts between the basal cell membrane and the cell substrate.

### System Includes:

- 16 and/or 96 well station located inside CO<sub>2</sub> incubator
- External control module
- Laptop PC
- ECIS control, acquisition, and display software
- Elevated field module (EFM) for automated cell migration and electroporation
- Twelve 8-well and/or six 96-well consumable electrode arrays




## 16 AND 96 WELL ARRAY STATIONS

ECIS® Array Stations connect to the back of the ECIS® Z and Zθ Station Controllers. They are usually placed inside a tissue culture incubator and have been designed to operate at high humidity conditions.

A flat shielded cable from the Array Station exits the incubator by the rubber seal of the inner glass door or through an incubator port, if available. LED's show the state of the device while looking through the incubator window.

The 16 Well Array Station provides electrical contact for two 8 well ECIS® arrays. The 96 Well Array Station provides electrical contact for a single 96 well ECIS® Array.

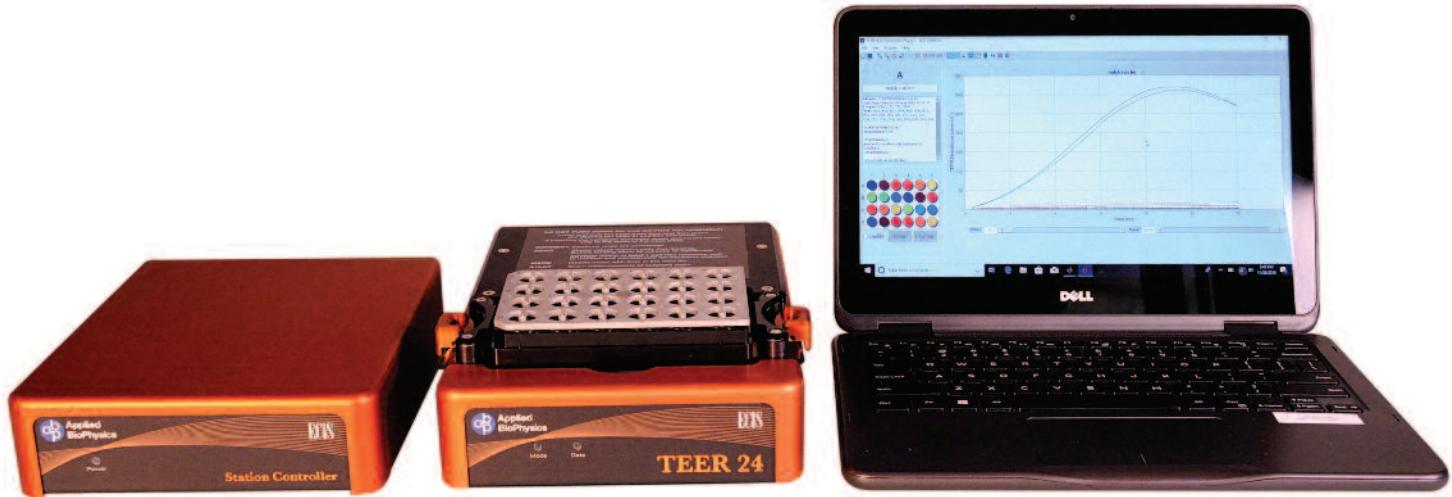

## TEER24

This system provides repeatable, label free automated TEER measurements to electrically monitor the barrier function of epithelial and endothelial cells as they are grown in normal CO<sub>2</sub> tissue culture incubators. Data is collected continuously and reported as real-time changes in barrier function of cell layers in 10-1000 ohm·cm<sup>2</sup>.

### System Includes:

- TEER 24 station located in CO<sub>2</sub> incubator
- Stainless steel plate assembly with 24 gold plated dipping pins
- External control module
- Laptop PC
- ECIS acquisition, control and display software
- Validation test array
- Four TEER24 consumable electrode arrays




## CP96 Cell Proliferation

The complete turn-key CP 96 system provides a means to carry out reproducible, label-free, automated cell proliferation measurements. Cell-proliferation is continuously monitored as cells grow in a normal CO<sub>2</sub> tissue culture incubator, and data are reported as real-time changes in percent cell coverage.


### System Includes:

- CP 96 station located in CO<sub>2</sub> incubator
- External control module
- Laptop PC
- ECIS acquisition, control and display software
- Validation test array
- Six 96-well consumable electrode arrays

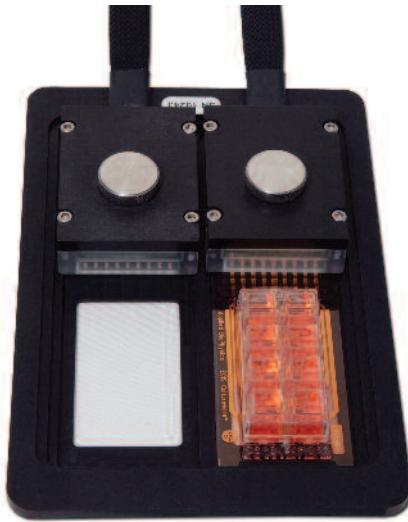

## ARRAY STAGE HOLDER

This device allows simultaneous ECIS® and optical measurements but is simply a platform and does not provide a chamber for heating and atmospheric control as with our Stage Incubator. The Array Stage Holder accommodates two 8 well ECIS® arrays (not compatible with ECIS® Flow Arrays) and fits on the stage of an inverted tissue culture microscope.

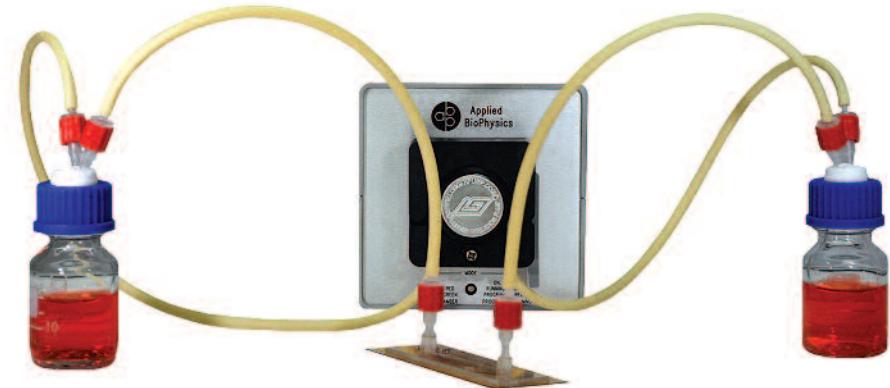

## P-FLOW PERISTALTIC PUMP

The Model ECIS p-Flow can be controlled manually to adjust flow rate from 0.20ml/min with a minimum flow rate of 70ul/min. Direction of the pump can be adjusted via a toggle switch or via the ECIS software (V 1.2.151 or higher).

The pump can also be controlled automatically from within the ECIS software (V 1.2.150 or higher) using USB or RS232. The pump can be programmed to run continuous or ramp up and down. It can store complex programmed flow profiles and run them disconnected from the PC. Each pump has a unique serial number so multiple pumps can be run from one PC.

For additional application information see Flow Based Assays, page 16.

### Performance specifications:

- Maximum Flow Rate with 3mm ID tubing: 20ml/min
- Minimum Flow Rate with 3mm ID tubing: 70ul/min (1.2uL/sec)
- Flow Rate resolution with 3mm ID tubing: 70ul/min from -20ml/min to +20ml/min
- Max ramp rate 0-20ml/min : 150msec (forward or reverse)
- Full forward to full reverse, -20ml/min to +20ml/min : 400msec
- Pump Speed: -60 to +60 RPM stepper driven motor
- Pump Resolution: 1024 increments/rev
- Control Input: USB or RS232



## CO<sub>2</sub> TISSUE CULTURE INCUBATOR

The ECIS® array holder is normally located in a tissue culture incubator to provide the normal growth requirements for cells. Since thermal and CO<sub>2</sub> fluctuation affect the behavior of the cells and can show up in ECIS® measurements, we recommend the use of a dedicated incubator that is only opened and closed occasionally during data acquisition.

The incubator has an exterior footprint of approximately 24 × 24 inches (28 inches in height) and a 1.8 cubic foot interior chamber. A port, with a special seal into the chamber, is provided for the ECIS® leads as well as for tubing if one is using the ECIS® flow module.

# 4 arrays

ECIS® Cultureware consists of patterned gold electrodes on thin plastic films bonded to different well configurations including 8 well arrays, 96 well microtiter plates, or specialized flow arrays from ibidi GmbH. The gold layer is thin enough to allow microscopic observation of the cells using a standard inverted tissue culture microscope



Applied  
BioPhysics

## STANDARD 96 WELL ARRAYS

ARRAY	ELECTRODES PER WELL	ELECTRODE AREA (mm <sup>2</sup> )	NUMBER OF CELLS MEASURED WHEN CONFLUENT	WELL VOLUME (μL)
96W1E+	2	0.256	100-200	300
96W10idf	idf	2.09	2000-40000	300
96W20idf	idf	3.985	4000-8000	300

## STANDARD 8 WELL ARRAYS

ARRAY	ELECTRODES PER WELL	ELECTRODE AREA (mm <sup>2</sup> )	NUMBER OF CELLS MEASURED WHEN CONFLUENT	WELL VOLUME (μL)
8W1E*	1	0.049	50-100	600
8W10E*	10	0.49	500-1000	600
8W10E+*	40	1.96	2000-4000	600
8WCP20idf	idf	3.985	4000-8000	600

## SPECIALTY ARRAYS

ARRAY	ELECTRODES PER WELL	ELECTRODE AREA (mm <sup>2</sup> )	NUMBER OF CELLS MEASURED WHEN CONFLUENT	WELL VOLUME (μL)
8W2x1E	2x1	2x0.049	50-100	600
8W1CXE	1	0.049	50-100	600
8W2LE	2	0.20	200-400	600
8Wμ1E+	4	0.196	200-400	600
2W4x10E	4x10	4x0.49	2000-4000	600

## FLOW ARRAYS

ARRAY	ELECTRODES PER WELL	ELECTRODE AREA (mm <sup>2</sup> )	NUMBER OF CELLS MEASURED WHEN CONFLUENT	CHANNEL/RESERVOIR VOLUME (μL)	CHANNEL HEIGHT X WIDTH (mm)
1F8x1E	8x1 (1 channel)	0.049	50-100	90/60	0.36 x 5
1F8x10E	8x10 (1 channel)	0.49	500-1000	90/60	0.36 x 5
6F1E	1 (6 channels)	0.049	50-100	45/60	0.66 x 5
6F10E	10 (6 channels)	0.49	500-1000	45/60	0.66 x 5
1F2Yx10E	8x4x2 (30&45 degree sides**)	0.49	500-1000	165/60	0.66 x 5

\*\*four electrodes are equal in size to 8W10E

\*available in PET and PC


idf: interdigitated fingers

PET: Polyethylene terephthalate, standard thickness 0.25mm

PC: Clear polycarbonate substrate, standard thickness 0.51mm or high numerical aperture (HNA) series of 0.13mm

PCB: Printed circuit board, standard thickness 1.55mm

## 96W STANDARD ARRAYS

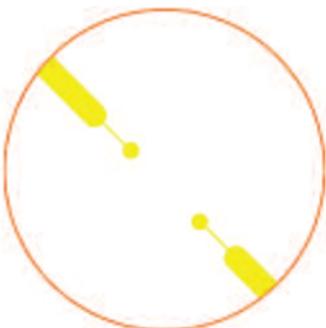


### 96W1E+

Each of the 96 wells in a standard plate configuration contains two circular  $350\mu\text{m}$  diameter active electrodes on a transparent PET substrate. As with other 1E arrays, a major use of this array is for the ECIS wound-healing assays where the small electrodes assure the high current pulse will result in complete cell killing.

Only a small population of cells is monitored on the small electrodes resulting in a fluctuating impedance signal due to the random like movement of the cells (micromotion).

#### Recommended Applications:

- Cell Migration
- In situ Cell Electroporation and Monitoring
- Measurement of micromotion
- Signal transduction assays

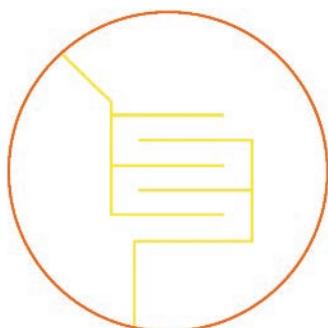

### 96W10idr

Each of the 96 wells has an Inter-digitated finger configuration. The total electrode area is  $1.96\text{mm}^2$  which measures 2000-4000 cells.

#### Recommended Applications:

- Signal transduction assays
- Detection of invasion of endothelial cell layers by metastatic cells
- Barrier function
- Cytotoxicity
- Cell differentiation




### 96W20idr

Each of the 96 wells has an Inter-digitated finger configuration. The total electrode area is  $3.92\text{mm}^2$  which measures a 4000-8000 cells.

#### Recommended Applications:

- Cell-ECM protein interactions
- Signal transduction assays
- Detection of invasion of endothelial cell layers by metastatic cells
- Barrier function
- Cell proliferation
- Cytotoxicity
- Cell differentiation

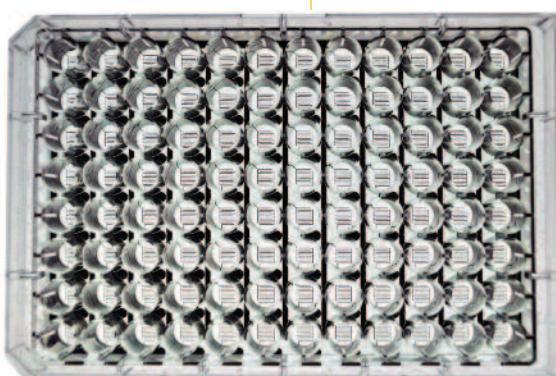

#### Array Color Key



Gold



Insulating Film

## 8W STANDARD ARRAYS



### 8W1E

Each well contains a single circular 250 µm diameter active electrode. On average, with a confluent cell layer, approximately 50 to 100 cells will be measured by the electrode, but even a single cell can be observed.

#### Recommended Applications:

- Signal Transduction
- In situ Cell Electroporation and Monitoring
- Cell Migration / Wound Healing
- Correlated microscopy and ECIS® experiments.




### 8W10E

Each well contains ten circular 250 µm diameter active electrodes connected in parallel on a common gold pad. On average, with a confluent cell layer, approximately 500 to 1000 cells will be measured by the electrodes.

#### Recommended Applications:

- Cell Differentiation
- Barrier Function
- Signal Transduction
- Cell Invasion
- Cytotoxicity

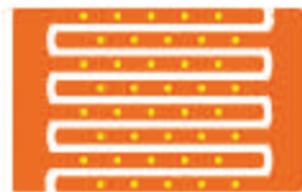

### 8W10E+

10E+: Each well has two sets of 20 circular 250 µm diameter active electrodes located on inter-digitated fingers to provide measurements of cells upon a total of 40 electrodes. On average, with a confluent layer, approximately 2000 to 4000 cells will be measured by the electrodes.


The 10E+ arrays are designed to monitor larger numbers of cells, sampling over the entire bottom of the well. Because of the relatively high number of cells, impedance fluctuations due to micromotion are largely smoothed out and do not obscure subtle changes in impedance due to the experimental conditions.



#### Recommended Applications:

- Cell Attachment and Spreading
- Cell Proliferation
- Cell Differentiation
- Cell-ECM Protein Interactions
- Barrier Function
- Signal Transduction
- Cell Invasion
- Cytotoxicity

## SPECIALTY ARRAYS

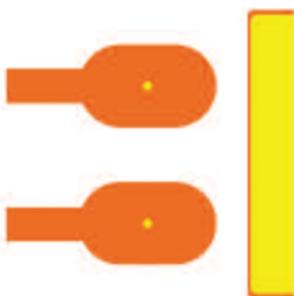

### 8W2x1E

This array is also called the Medusa array. Each well in this array has two independent single 250 µm diameter active electrodes. The Medusa array is useful for duplicating readings in the same well or to wound/electroporate one electrode while leaving the other as a control within the same well.

When connected to the array holder only the upper four wells are measured. To use the other four wells, the array is turned around and the contact pads at the other end are connected.

#### Recommended Applications:

- Signal Transduction
- In situ Cell electroporation and Monitoring
- Cell Migration / Wound Healing
- Correlated microscopy and ECIS® experiments.




### 8W20idf

Our special purpose cell proliferation array samples the bottom of the 8 well chamber with our patterned electrodes. As few as 5000 cells per well can be seeded and detected by this array allowing for cell proliferation to be monitored over the course of approximately 5 cell doubling times.

The 8WCP arrays are designed to monitor larger numbers of cells, sampling over the entire bottom of the well. Because of the relatively high number of cells, impedance fluctuations due to micromotion are almost completely smoothed out and do not obscure subtle changes in impedance due to the experimental conditions.

#### Recommended Applications:

- Cell Attachment and Spreading
- Cell Proliferation
- Cytotoxicity

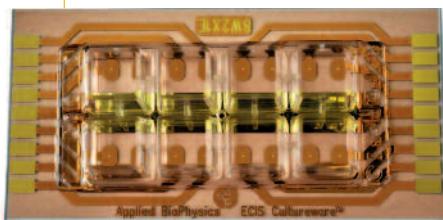

#### Array Color Key

Gold

Insulating Film



### 8W1CXE

Also known as the Chemotaxis array, its use was first described in 2001\*. The measuring electrode in this array is a thin gold line 75 µm x 667 µm between two registry marks. The gold line has the same total area as a 250 µm single circular electrode.


In the ECISaxis assay described by Hadjout et al. the well is flooded with agarose and allowed to polymerize. Once the gel has hardened, wells in the agarose are introduced above the circular registry marks. Chemoattractant is then added in one well, while cells are added to the other. The single electrode is sensitive enough to detect the migration of single cells.

\*by Hadjout, N. et al. (2001)  
*Biotechniques* 31 (5) 1130.

#### Recommended Applications:

- Cell Chemotaxis

## SPECIALTY ARRAYS

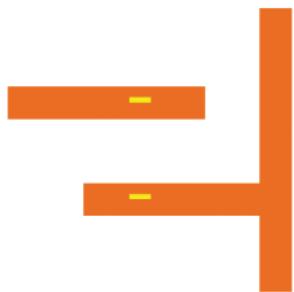


### 8W2LE

Each of the 8 wells contains two linear electrodes with dimensions of 667 $\mu\text{m}$  x 150 $\mu\text{m}$  in series. Each well thus has an electrode area four times that of our standard 250 $\mu\text{m}$  circular electrodes but retains the same impedance values. These arrays were designed for cell migration measurements in which rectangular cell-free areas are generated for direct comparisons with traditional scratch assays.

#### Recommended Applications:

- Cell Migration / Wound Healing
- Correlated microscopy and ECIS® experiments

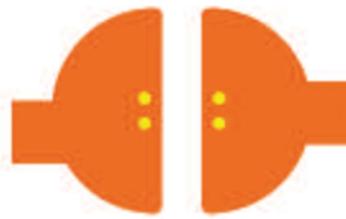

### 8W $\mu$ 1E+

Each of the 8 wells contains four 250 $\mu\text{m}$  circular electrodes. The placement of the electrodes at the center of the well allows for the use of cloning cylinders to be placed around the electrodes creating micro-wells. The area outside of the cloning cylinder can then be flooded to reduce evaporation from within the micro-wells.

#### Recommended Applications:

- Barrier
- Signal Transduction
- Cell Invasion
- In situ Cell electroporation and Monitoring
- Cell Migration / Wound Healing
- Correlated microscopy and ECIS® experiments




### 2W4x10E

Each of the 2 circular 25 mm diameter wells contains four independent sets of ten 250  $\mu\text{m}$  diameter active electrodes. This array was specifically designed for hypoxia studies to create a large liquid – air interface for rapid gas exchange. It's design also incorporates a small central area devoid of gold or photoresist allowing for live cell fluorescence microscopy.

#### Recommended Applications:

- Cell Attachment and Spreading
- Cell Proliferation
- Cell Differentiation
- Barrier Function
- Signal Transduction
- Cell Invasion
- Cytotoxicity
- Correlated microscopy and ECIS® experiments

## FLOW ARRAYS

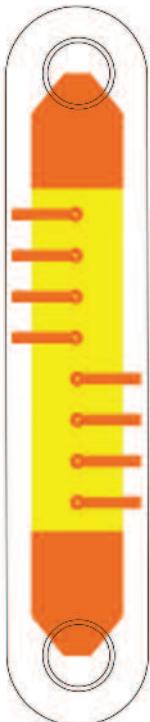


### 1F8x1E and 1F8x10E

These are specialized Flow arrays with the 1F8x1E having 8 active 250 µm diameter electrodes, and the 1F8x10E having 8 sets of 10 active 250 µm diameter electrodes located in the central region at the base of a flow channel measuring 50 mm in length 5 mm in width and available in 0.36 mm in height with a total channel volume of 90 µL.

Our flow arrays are designed for ECIS® measurements of cells under perfused conditions or under flow mimicking the shear stress endothelial cells experience *in vivo*.

#### Recommended for the following applications under shear stress conditions:

- Barrier Function
- Signal Transduction
- Cell Invasion
- *In situ* Cell electroporation and Monitoring
- Cell Migration / Wound Healing
- Cell Proliferation
- Cell Differentiation
- Barrier Function
- Cytotoxicity

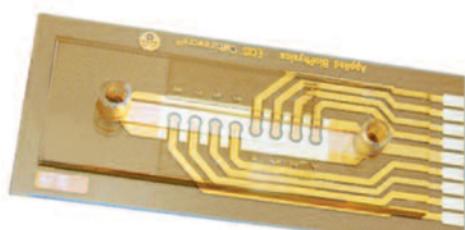

Array Color Key





Gold

Luer lock connection



Insulating Film

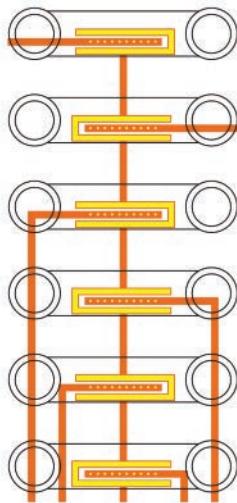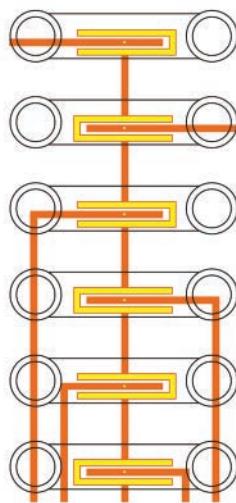


### 6F1E and 6F10E PC


This flow array allows 6 independent flow assays to be run simultaneously. The channels are 0.66mm in height and 5mm wide with either 1 or 10 active 250µm diameter electrodes per channel. Each channel has a 45µL volume with 60µL reservoirs.

#### Recommended for the following applications under shear stress conditions:

- Barrier Function
- Signal Transduction
- Cell Invasion
- *In situ* Cell electroporation and Monitoring
- Cell Migration / Wound Healing
- Cell Proliferation
- Cell Differentiation
- Barrier Function
- Cytotoxicity

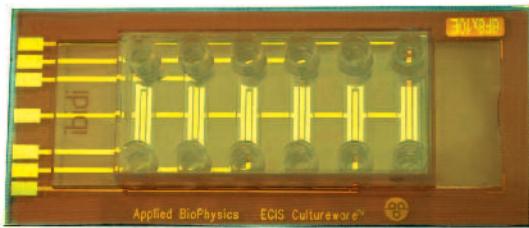


Luer lock connection

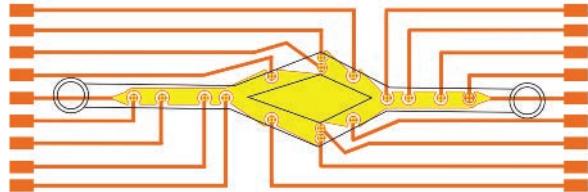

### 1F2Yx10E PC


This flow array is intended for bifurcation studies and blood vessel simulation. It splits into 30 degree Y channels in one direction and 45 degree Y channel in the other direction.

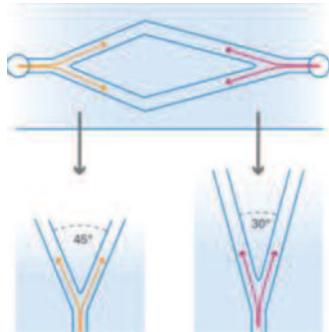

This array is double ended with 8 measurement channels available at each end. Eight measurement points, each with 4 circular active electrodes (which are equal in area to an 8W10E), are located along the channel and through the Y portion of the channel. One end of the array is used to monitor the 30 degree Y channel and the other end is used to monitor the 45 degree Y channel. The electrodes are located close in the corners of the flow direction transition points. Each channel has a 165 $\mu$ L volume with 60 $\mu$ L reservoirs. The flow is always laminar, i.e., turbulent flows are not possible. For simulation of turbulence flow we recommend oscillating the flow. Defined shear stress and shear rate levels.

#### Recommended for the following applications under shear stress conditions:

- Simulation of the bifurcation of blood vessels for arteriosclerosis research
- Rolling and adhesion of leukocytes on endothelial cells cultured under flow
- Cell-cell interaction studies and cell-drug interaction screenings under flow conditions

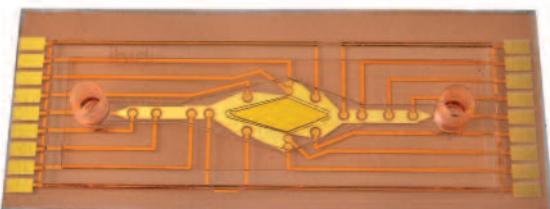

Luer lock connection

### CUSTOM ARRAYS

If you are interested in other electrode sizes and configurations please contact us. We will be happy to work with you on your specific needs.

### HOW TO ORDER ARRAYS

We maintain a supply of arrays for shipment from our facility. The arrays are shipped sterile.


All arrays can be ordered without upper chambers (NC). Gold films are transparent and the upper chambers can be removed following ECIS® measurements to facilitate microscopy including fluorescence measurements

We recommend arrays not to be stored for more than two months.

For ordering, please contact 518-880-6860 or send email to [info@biophysics.com](mailto:info@biophysics.com).

Please include shipping and billing address.

Credit cards are accepted.

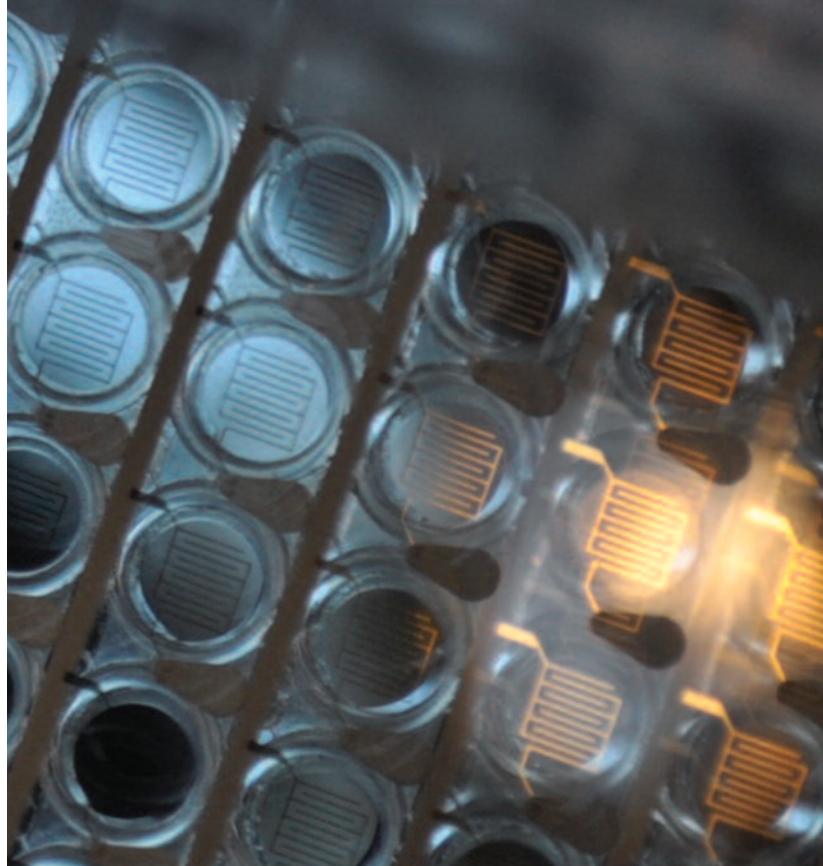

# institutions that use ECIS

## Australia

University of Sydney  
University of South Australia

## Austria

Biotec Area Krems  
Medical University Graz  
Tissue Med Biosciences GmbH & Co  
University of Applied Sciences Krems

## Belgium

Ghent University  
Katholieke Universiteit Leuven

## Brazil

Natura Invacao e Technological

## Canada

BioPhage Inc  
Conseil National de Recherches  
Fisheries and Oceans Canada  
Millenium Biologix  
Ottawa Hospital  
St. Michael's Hospital  
Université Laval  
University of Alberta  
University of Montreal  
Vertex Pharmaceuticals Incorporated

## China

Capital University of Medical Sciences  
Central Southern University  
Fu Jen Catholic University  
Fudan University Medical School  
Guangdong Entry-Exit Inspection and Quarantine Bureau  
Guangzhou Medical University (GMU)  
Inner Mongolia Medical College  
Institute of Materia Medica  
Jinan University  
Shandong Provincial Qianfoshan Hospital  
Shanghai Jiao Tong University  
Sun Yat-sen University  
The Third Xiangya Hospital of Central South University  
ZhenJiang First People's Hospital

## Denmark

Novo Nordisk

## Estonia

Tallinn University of Technology

## Finland

Åbo Akademi University  
University of Oulu

## France

Centre Hospitalier Universitaire de Nice, Inserm, U576  
Centre de recherche Cardiovasculaire à l'HEGP, Inserm, U970

## Germany

Albert-Ludwigs Universität  
Bayer Schering Pharma AG  
Beiersdorf AG  
Charite Universitätsmedizin Berlin  
CLR Chemisches Lab.  
Fraunhofer Institute for Interfacial Engineering & Biotechnology IGB  
Friedrich-Alexander-Universität  
Friedrich-Schiller-University Jena  
Georg-August-Universität Göttingen  
Johannes Gutenberg University Mainz  
Johann Wolfgang Goethe-University  
Ludwig-Maximilians-Universität München  
Max-Delbrück-Centrum für Molekulare Medizin (MDC)  
MetaGen Pharmaceuticals GmbH  
Philipps-Universität Marburg  
Saarland University  
Technical University Munchen  
Universitat Regensburg  
Universitat Wurzburg  
University Medical Centre Mannheim  
University of Tuebingen  
University of Witten/Herdecke  
WWU Munster

## Hungary

Debreceni Egyetem  
Semmelweis University  
University of Debrecen Medical School

## India

Bengal Engineering and Science University, Shibpur  
Central Leather Research Institute  
Indian Institute of Science Bangalore  
Indian Institute of Technology Madras  
Jawaharlal Nehru Centre for Advanced Scientific Research (JNCASR)  
National Centre for Cell Science

National Institute of Immunology  
Osmania University

## Ireland

Trinity College Dublin  
University College Cork

## Israel

Ben Gurion University of the Negev  
Haifa University  
Rambam Hospital  
Teva Pharmaceutical Industries Ltd.

## Italy

Analitica DeMori Srl.  
European Commission Joint Research Centre  
IFOM-IEO Campus

## Japan

Astellas Pharma Inc.  
Fukui Medical School of Japan  
Kawasaki Medical School  
National Cancer Center  
National Center for Child Health and Development  
National Center for Global Health and Medicine Research Institute  
National Cerebral and Cardiovascular Center  
Okayama University  
Osaka City University Medical School  
Otsuka Pharmaceutical Co.  
University of Tokyo  
University of Tsukuba

## Korea

DaeJeon Science High School  
Chungbuk National University  
Gachon University  
Jeju National University  
Kyungwon University  
Yonsei University College of Medicine

## Malawi

Queen Elizabeth Central Hospital

## Malaysia

University of Malaya

## Netherlands

AMC University of Amsterdam  
Leiden University Medical Center  
Sanquin Research CLB

## UMC Utrecht

University Medical Center Groningen  
VU University Medical Center

## Norway

University of Bergen  
University of Oslo

## Qatar

Qatar University

## Russia

Russian Cardiology Research Center

## Spain

Universitat Autonoma de Madrid

## Sweden

AstraZeneca  
Karolinska University Hospital, Solna  
Umeå University

## Switzerland

Merck Serono  
Novartis Pharma AG  
University Hospital Zurich

## Taiwan

Academia Sinica Institution  
Chang Gung University  
Chi Mei Hospital  
Food Industry R and D Institute  
Fu Jen Catholic University  
National Cheng Kung University  
National Defense Medical Center  
National Chiao Tung University  
National Taiwan University  
National Taiwan University Hospital  
National Yang-Ming University  
Tzu Chi College of Technology

## United Kingdom

Cardiff University  
Novartis Institute for BioMedical Research  
Queen's University  
University College London  
University of Bristol  
University of Edinburgh  
University of Manchester  
University of Portsmouth

**United States**

Abbott Laboratories  
 Al Dupont Hospital for Children  
 Albany Medical College  
 Albert Einstein College of Med.  
 Alcon Laboratories  
 Allergan  
 Amgen Inc.  
 Armed Forces Radiobiol. Res. Inst.  
 Bausch & Lomb  
 Baylor College of Medicine  
 Beth Israel Deaconess Medical Center  
 Blood Center of Wisconsin  
 Blood Systems Research Ins.  
 Boston Children's Hospital  
 Brigham and Women's Hospital  
 Bristol-Myers Squibb  
 Buffalo State  
 Capstone Therapeutics Corp.  
 Case Western Reserve University  
 Celgene Corporation  
 Cellular Dynamics International Inc.  
 Centers for Disease Control and Prevention  
 Charlie Norwood VA Medical Center  
 City College of New York  
 Cleveland Clinic Foundation  
 Columbia University  
 Department of Veterans Affairs Medical Center  
 Eli Lilly and Company  
 Emory University  
 Florida International University  
 Georgetown University  
 Georgia Regents University  
 Georgia State University  
 HemoShear LLC  
 Henry Ford Hospital  
 Howard Hughes Medical Center  
 Indiana University  
 IUPU  
 Indianapolis  
 Johns Hopkins University  
 Johnson & Johnson  
 Kean University  
 La Sierra University  
 LSU School of Medicine  
 Mary Kay Inc.  
 Mayo Clinic Jacksonville  
 Mayo Clinic Rochester  
 MD Anderson Cancer Center  
 Medical University of South Carolina  
 MedImmune Inc.

Merck & Co.  
 Minnesota State University Moorhead  
 Morehouse School of Medicine  
 National Cancer Institute  
 National Institute of Environmental Health Science (NIEHS)  
 National Institutes of Health (NIH)  
 National Institute for Occupational Safety and Health (NIOSH)  
 North Carolina Agricultural & Technical State University  
 Northwestern University  
 Novo Nordisk  
 NYU Langone Medical Center  
 Ohio State University  
 Omaha VA Medical Center  
 Oregon Health & Science University  
 Oregon Med Laser Center  
 Procter & Gamble Co.  
 Providence VA Medical Center  
 Rensselaer Polytechnic Institute  
 Rhode Island Hospital  
 Roswell Park Cancer Institute  
 Rush University Medical Center  
 Seattle Children's Hospital  
 Spelman College  
 St. Louis University  
 St. Michael's Hospital  
 Stanford University  
 Stemnion, Inc.  
 SUNY Cortland  
 Temple University School of Medicine  
 Texas A&M University  
 The Children's Mercy Hospital  
 The Commonwealth Medical College  
 Trauma Research LLC  
 Trinity Biosystems  
 U.S. Food and Drug Administration  
 University of Alabama  
 University of California Davis  
 University of California Irvine  
 University of California Los Angeles  
 University of California San Francisco  
 University of Chicago  
 University of Cincinnati  
 University of Connecticut  
 University of Delaware  
 University of Florida Gainesville  
 University of Hawaii at Manoa  
 University of Illinois at Chicago  
 University of Kentucky  
 University of Louisville  
 University of Michigan

University of Minnesota  
 University of Mississippi Medical Center  
 University of New Mexico  
 University of Oklahoma Health Sciences Center  
 University of Pennsylvania  
 University of Pittsburgh  
 University of Rhode Island  
 University of Richmond  
 University of South Alabama College of Medicine  
 University of South Carolina  
 University of South Florida Tampa  
 University of Tennessee Knoxville  
 University of Texas  
 University of Texas Health Science Center at San Antonio  
 University of Texas Medical Branch  
 University of Texas Southwestern Medical Branch  
 University of Utah  
 University of Vermont  
 University of Virginia  
 University of Wisconsin  
 US Army Center for Environmental Health Research  
 US Army Dental & Trauma Research Detachment  
 US Army Institute of Surgical Research  
 US Army Medical Research Institute of Infectious Disease  
 VA Greater Los Angeles Healthcare System  
 VA Medical Center, Dallas  
 Van Andel Research Institute  
 Vanderbilt University  
 Virginia Commonwealth University  
 Walter Reed Army Institute of Research  
 Wayne State University  
 Weill Cornell Medical College  
 West Virginia University  
 Womens & Infants Hospital  
 Yale University

# international distributors

**ECIS®** is distributed world wide with instruments on every continent except Antarctica. Please contact us through the distributor in your region.

## AUSTRALIA

**SDR Clinical Technology**  
206 / 354 Eastern Valley Way  
Chatswood NSW 2067  
Tel: 61-2-9882-2882  
Fax: 61-2-9882-6468  
Attn: Dr. Peter Kenny  
Web: <http://www.sdr.com.au/>

## CHINA

**Dakewe Biotech Co**  
Unit B, 3/F, The 5th Building  
Tian'an industrial village  
Nanyou, Dengliang Rd., Nanshan Dist.,  
Shenzhen 518054 China  
Tel: 0086 755 2686 946  
Fax: 0086 755 2686 9149  
Attn: Dr. Qing Jun  
Web: [www.dakewe.net](http://www.dakewe.net)

## ESTONIA

**Inbio OÜ**  
Akadeemia tee 15  
Tallinn 12618  
Tel: +372 6204350  
Fax: +372 6204358  
Attn: Harri Peiker  
Email: [harri@inbio.ee](mailto:harri@inbio.ee)  
Web: [inbio.ee](http://inbio.ee)

## GERMANY

**ibidi GmbH**  
Am Klopferspitz 19  
D-82152 Martinsreid (Munchen)  
Tel: 49 89 5204 61731  
Fax: 49 89 5204 61759  
Attn: Dr. Ulf Radler  
Email: [uraedler@ibidi.de](mailto:uraedler@ibidi.de)  
Web: <http://www.ibidi.de/>

## INDIA

**Labmate**  
Biad Mehta Complex  
183 Mount Road  
Chennai, India 600 015  
Tel: 44 2220 0066  
Fax: 44 2220 0056  
Attn: S. Viswanathan  
Email: [labmate@vsnl.com](mailto:labmate@vsnl.com)  
Web: <http://www.labmateasia.com/>

## ISRAEL

**Almog Diagnostic**  
PO Box 634  
Rishon Le Zion  
Israel 75103  
Tel: 972 3 967 3390  
Fax: 972 3 967 3391  
Attn: Nitsan Levi  
Email: [nitsan@almog.co.il](mailto:nitsan@almog.co.il)  
Web: [www.almog.co.il](http://www.almog.co.il)

## ITALY

**SisLab Sas**  
Via Martiri Oscuri 10  
20125 Milan – Italy  
Tel: 39 02 222 23968  
Fax: 39 02 928 76689  
Attn: Riccardo Sacco  
Email: [info@sislab.it](mailto:info@sislab.it)  
Web: <http://www.sislab.info>

## JAPAN

**Nepa Gene Company, Ltd.**  
3-1-6 Shioyaki  
Ichikawa, Chiba 272-0114 Japan  
Tel: 81 47 306 7222  
Fax: 81 47 306 7333  
Attn: Suzuki  
Web: [www.nepagene.jp](http://www.nepagene.jp)

## KOREA

**LeeBaeg Science Co., Ltd.**  
Second Floor  
626-1 Guam-Dong  
YuSung-Gu, Daejeon-City 305-801  
Email: [info@lpscience.com](mailto:info@lpscience.com)  
Web: <http://www.lpscience.com>

## MALAYSIA

**Precision Technologies**  
211 Henderson Road #13 - 02  
Henderson Industrial Park  
Singapore 159552  
Tel: (65) 6273-4573  
Fax: (65) 6273-8898  
Attn: Lim Kian Heng  
Web: <http://www.pretech.com.sg>

## PUERTO RICO

**Audio & Video Pro Caribe**  
J.C. Gonzalez, Inc.  
HC 01 Box 8389  
San German, P.R. 00683-9714  
Tel: 787-892-0047  
Fax: 787-264-3816  
Attn: Julio C. Gonzalez  
Email: [jcgonzalez@avpcaribe.com](mailto:jcgonzalez@avpcaribe.com)  
Web: [www.avpcaribe.com](http://www.avpcaribe.com)

## SINGAPORE

**Precision Technologies**  
211 Henderson Road #13 - 02  
Henderson Industrial Park  
Singapore 159552  
Tel: (65) 6273-4573  
Fax: (65) 6273-8898  
Attn: Lim Kian Heng  
Web: <http://www.pretech.com.sg>

## TAIWAN


**Sunpoint Scientific Instrument Co., Ltd.**  
9F-a, No. 146, Wen Shing Road  
Kwei Shan Hsiang  
Taoyuan Hsien 333  
Tel: 886 33 27388  
Attn: Sunny Tsai  
Email: [sunp@ms36.hinet.net](mailto:sunp@ms36.hinet.net)  
Web: <http://www.sunpointworld.com/>

## TURKEY

**aka Biyoteknoloji**  
Nazmi Akbacı Ticaret Merkezi  
No.246 Maslak 34398 ISTANBUL  
Tel: +90(0)212 346 2019  
Attn: Alp Kibaroglu  
Email: [alpkibaroglu@akabiotech.com](mailto:alpkibaroglu@akabiotech.com)  
Web: [www.akabiotech.com](http://www.akabiotech.com)

## UNITED KINGDOM

**ibidi GmbH**  
Am Klopferspitz 19  
D-82152 Martinsreid (Munchen)  
Tel: 49 89 5204 61731  
Fax: 49 89 5204 61759  
Attn: Dr. Ulf Radler  
Email: [uraedler@ibidi.de](mailto:uraedler@ibidi.de)  
Web: <http://www.ibidi.de/>

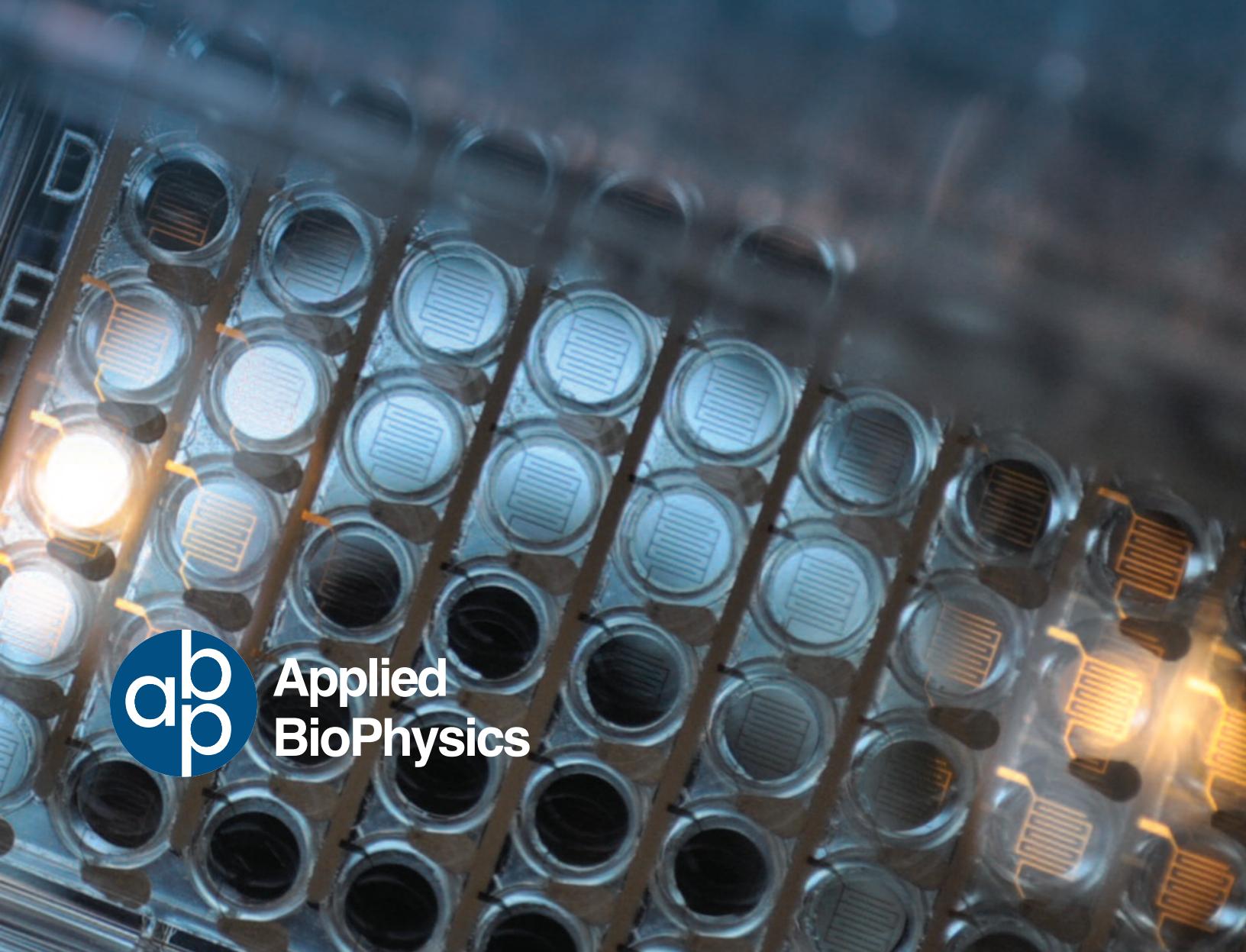

Applied  
BioPhysics

## Request A Demo:

Applied Biophysics (ABP) agrees to make available an ECIS instrument for a 4-week period.

All equipment necessary for Electric Cell-substrate Impedance Sensing is provided by ABP, which includes: ECIS controller, array station, computer and start up supply of arrays. Additional arrays are available for purchase.

ABP agrees to pay the shipping from ABP to your location and should the decision be made not to keep the instrument, requests the customer pays return shipping fees.

ABP requests at least one researcher be dedicated to the operation and maintenance of the instrument. As the dedicated operator we request to speak with this researcher on a weekly basis to review data and system operation.

To request a demo please contact [info@biophysics.com](mailto:info@biophysics.com) or contact customer support @ **1-866-301-3247**

## Webinars:

The ECIS Application Webinar series will review the topics listed below in 20-minute, web-based, interactive seminars.

- ECIS Theory
- Cell Invasion/Extravasation Assays
- Automated Wounding and Cell Migration
- Barrier Function Assays
- Real-time Electroporation and Monitoring
- Cell Attachment and Spreading Measurements
- Signal Transduction Assays
- Toxicology with ECIS

For more information visit: [www.biophysics.com/webinar.php](http://www.biophysics.com/webinar.php)

## ECIS School:

Applied BioPhysics offers a 2-day in-depth training class several times per year.

The course, taught by Dr. Charles Keese, provides the students with the knowledge and tools on how to apply ECIS theory in cellular research and how to properly operate the ECIS systems to further enhance potential publication opportunities.

Applied BioPhysics will also provide an opportunity to get familiar with our area by arranging lunch and dinner at local eateries for the duration of your stay.

For more information visit: [www.biophysics.com/school.php](http://www.biophysics.com/school.php)





# Applied BioPhysics

Corporate Headquarters:  
185 Jordan Road • Troy, NY 12180  
1-866-301-ECIS (3247)

[www.biophysics.com](http://www.biophysics.com)

